
WELCOME

Congratulations on your new Hypervolt

You are one step closer to moving better, feeling better, and being a better version of you. The Hypervolt was created to keep you doing what you love. Whether you're an elite athlete, weekend warrior, or simply trying to get your body to feel its best, the Hypervolt will help you get there... and then some.

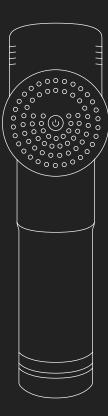

NEXT STEPS

Let's get started

1

Register your device at hyperice.com/register-product to activate your warranty and ensure easy returns, repairs, or refunds.

2

Download the Hyperice App to get the most out of the Hypervolt, or any Hyperice device.

App Store

3

Charge your device by plugging in your new Hypervolt with the provided Hyperice charger. Fully charge the battery for up to six hours before first use.

Have questions? Call 1-855-734-7224, chat live at hyperice.com, or email customersupport@hyperice.com.

QUICK START GUIDE

Moving better, made simple



Power up

Switch your Hypervolt to the ON position at the base of the handle. Once your device is powered ON a blue light will illuminate, plus (up to) four green lights for battery level.





Activate pressure sensor


Activate the pressure sensor by pressing the speed setting button once. Wait five seconds as the white lights sequence and turn off. **Your Hypervolt is now ready for use.**





Adjust speeds


Choose your desired speed (1-3) by pressing the speed setting button once per speed. The blue lights above the speed setting button will indicate the current speed.

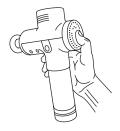

5

Recharging

Only use the provided Hyperice wall charger to uphold warranty. Full charge is indicated when the LED "band" on the handle stops flashing and becomes solid green.

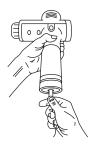



Changing attachments


Choose desired head attachment and insert directly into the shaft / opening while pushing firmly. **Fork attachment?** Align the "groove" on the fork with the recess on the shaft / opening and push firmly.






Caring for your Hypervolt


Use a damp, clean, cloth and gently wipe your Hypervolt. Make sure the power is OFF and the battery charger is not attached. Store in a clean, cool place when not in use.

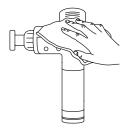

Choosing the best head attachment

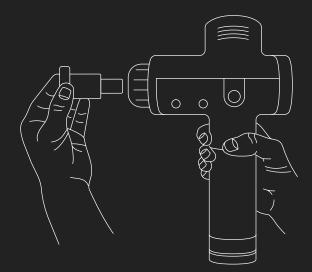

HEAD ATTACHMENTS

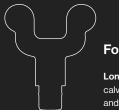

Fork

Longer muscles such as calves, forearms, biceps, and ankle / Achilles area.



Trigger points found in glutes, hips, traps, shoulders, and other complex muscle groups.



Flat

Larger surface areas such as quads, hamstrings, back, and chest.



Cushion

Sensitive areas such as wrists, ankles, knees, traps, and around the neck.



Ball

Sore muscles and a great starter attachment for larger muscle groups.

Want to learn more? Visit hyperice.com/how-to for video tutorials on your Hypervolt Plus or any Hyperice device.

HYPERICE APP

Your digitally connected wellness guide

Connect, monitor progress, learn, and operate Hyperice products from your IOS devices.



HYPERICE APP

A few of the features



Apple Health[™] and Strava[™] connected

Connect to like-minded apps for a more personalized digital experience.



Data tracking / analytics

Tap into activity-based analytics to monitor your progress.



Personalized programs

Guided warmup, recovery, and maintenance sessions.



Operate your devices

Operate your Hypervolt or any Hyperice Bluetooth^{*} connected device.

🛞 Bluetooth



IOS is a registered trademark of Cisco and is used under license. The Apple logo, Apple Health, and the App Store are registered trademarks of Apple, Inc: and used under license. Strava is a registered trademark of Strava, Inc. and is used under license. The Bluetooth® word mark and logos are a registered trademark owned by Bluetooth SIG, Inc. and is used under license.

HYPERCARE

A dedicated team just for you

Every Hyperice product is backed by a 1-year warranty and award-winning support from our Hypercare[™] team – a group of specialists dedicated to your overall well-being and expert guidance on Hyperice products.



HYPERCARE

Here's a few of the things your team can do for you



Product knowledge

Answer questions about your device and best practices.

	())	
C	_]
1				Γ
		- 1		L

Place new orders

Help with placing orders online and/or authorized resellers.



Track new orders

Track new orders and work with carriers to find solutions.

ſ	\$

Billing solutions

Research and provide solutions for any billing issues.



Trouble shoot

Help trouble shoot and solve issues with your device.



Initiate returns

Provide instructions for returns, refunds, or repairs.


Get in touch.

Hypercare® is registered trademark of Hyperice, Inc. Copyright 2020 © Hyperice, Inc. All rights reserved. Call 1-855-734-7224 Chat hyperice.com Email customersupport@hyperice.com Hyperice, Inc 525 Technology Drive, Suite 100 Irvine, California 92618


A family of products designed to upgrade the human body

Hypervolt Plus

Percussion Massage Device

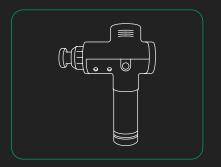


Vyper 2.0 Vibrating Fitness Roller

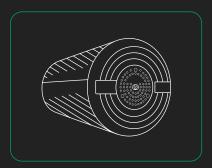


Hypersphere

Vibrating Massage Ball

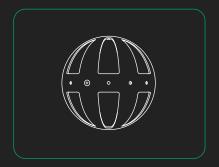

Venom

Wearable Heat + Vibration Devices

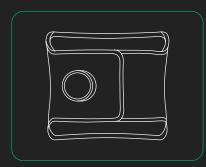


Hypersphere Mini Compact Vibrating Massage Ball



ICT

Ice Compression Devices

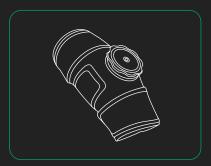

Learn more / shop @ hyperice.com

Free Manuals Download Website <u>http://myh66.com</u> <u>http://usermanuals.us</u> <u>http://www.somanuals.com</u> <u>http://www.4manuals.cc</u> <u>http://www.4manuals.cc</u> <u>http://www.4manuals.cc</u> <u>http://www.4manuals.com</u> <u>http://www.404manual.com</u> <u>http://www.luxmanual.com</u> <u>http://aubethermostatmanual.com</u> Golf course search by state

http://golfingnear.com Email search by domain

http://emailbydomain.com Auto manuals search

http://auto.somanuals.com TV manuals search

http://tv.somanuals.com