DK	Wellness lampe	2
SE	Harmonilampa.....	4
NO	Lysterapilampe.....	6
FI	Kirkasvalolamppu.....	8
UK	Wellness Lampe	10
DE	Wellnesslampe	12
PL	Lampa lecznicza do terapii światłem	14

INTRODUKTION

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INTRODUCTION

To get the best out of your new appliance, please read this user guide carefully before using it for the first time. We also recommend that you keep the instructions for future reference, so that you can remind yourself of the functions of your appliance.

SAFETY MEASURES

Normal use of the appliance

- Incorrect use of this appliance may cause personal injury and damage.
- Use for its intended purpose only. The manufacturer is not responsible for any injury or damage resulting from incorrect use or handling (see also Guarantee Terms).
- For connection to 230 V, 50 Hz only.
- Do not submerge in water or any other liquid.
- Do not use in humid surroundings.
- Do not touch with wet hands.
- Never leave unattended when in use, and keep an eye on children.
- Not suitable for commercial or outdoor use.
- Do not cover.

Cord, plug and mains socket

- Always turn the appliance off using the on/off button, and unplug it, when it is not in use.
- Keep cord away from hot objects and open flame.
- Check regularly that neither the cord nor plug is damaged and do not use if there is any damage, or if it has been dropped on the floor, dropped in water or damaged in any other way.
- If the apparatus or plug is damaged, it must be inspected and if necessary, repaired by an authorised repair engineer, otherwise there is a risk of electric shock. Never try to repair the appliance yourself.
- Remove the plug from the socket before cleaning.
- Avoid pulling the cord when removing the plug from the socket. Instead, hold the plug.
- Do not wind the cord around the appliance.
- Check that it is not possible to pull or trip over the cord or any extension cord.

Use of light therapy

- Do not use in rooms without some other light source. For use as a light supplement only, not illumination.
- Do not look directly at the apparatus for long periods. Follow the user guide below.
- Initial use can cause tired eyes and/or headache. This is normal, and such effects will soon disappear. Seek medical advice if they do not.
- Consult your doctor before using this apparatus if you
 - o suffer from or have suffered from severe depression
 - o are over-sensitive to light
 - o have painful eyes
 - o use particular forms of medicine such as psychopharmacia or malaria pills
 - o suffer from high blood pressure.
- Consult your doctor if the apparatus causes physical or mental symptoms of any sort, or if you have any other medical questions or problems.

KEY TO MAIN COMPONENTS OF THE APPLIANCE

1. Screen
2. Wall-hanging hole
3. On/off button
4. Carrying handle
5. Screen retaining screws
6. Stand
7. Safety bracket
8. Light tubes
9. Starter
10. Socket

USING THE APPARATUS

For light baths/light therapy to counter the effects of winter blues, which can cause the following symptoms:

- Tiredness
- Bad moods or mood swings
- lack of energy
- Increased need for sleep
- Inability to concentrate
- Loss of appetite

Over 60% of the population are believed to experience one or more symptoms of winter blues due to the low intensity of the winter sun.

The most effective time of day to take a light bath is between 6 am and 8 pm for a recommended period of 2 hours per day. You can take a light bath as often as you like, but treatment is most effective when you take one for the times below for at least 5 successive days.

Positioning the appliance

- Fold the stand (6) out and place the apparatus on a table top. or
- Hang the apparatus on a wall using the hole (2) on the back.
- o Check that the wall and screw/nail can bear the weight of the apparatus before hanging it up.

Use

- Connect the plug to the mains socket and switch on using the on/off switch (3).
- Sit in front of the apparatus so that you can see the light easily and with your face 20-60 cm from the screen (1). You do not need to sit in front of the screen.
- How long your light bath should last depends on the distance you sit from the screen:

Distance to screen	Max. length of light bath
60 cm	Approx. 2 hours
40 cm	Approx. 1 hour
20 cm	Approx. ½ an hour

- You can take a break during the session and do other things whilst bathing (e.g. read or watch TV) as long as you sit within 60 cm of the apparatus.
- Do not look directly into the light constantly, but look at it regularly.
- Always turn the appliance off using the on/off button, and unplug it, when it is not in use.

Allow to cool for at least 15 minutes before putting it away or carrying it by the handle (4).

REPLACING LIGHT TUBES

Normal tube (8) lifetime is approx. 10,000 hours of normal use. If a tube fails, it must be replaced.

- Slacken the retaining screws (5) on the back using a screwdriver, lay the apparatus on its back and remove the screen.
- o Note that the four brackets on the front of the screen will become loose when the screws are slackened. Take care not to lose them, as the screen cannot be refitted correctly.
- Lift the safety bracket (7) and carefully lift upper end of tube from the apparatus.
- Pull the tube free of its socket (10) at the bottom of the apparatus. Dispose of used tubes in accordance with local rules.
- Insert a new tube in the reverse order.
- o Both tubes should be replaced at the same time.
- Replace the screen.

CLEANING

- Remove the plug from the socket before cleaning.
- Allow to cool down before cleaning.
- Do not use strong or abrasive cleaning agents, solvents or similar to clean the apparatus, as these may scratch or damage the external surfaces.
- Never submerge the apparatus in water, and no water must be allowed to get into it.
- Allow to dry fully after cleaning before using again.

External cleaning

- Wipe off with a damp, soft cloth.

Internal cleaning

- Remove the screen as described above, wipe off reflector, inside of screen and tubes with a damp, soft cloth. Replace the screen.

INFORMATION ON DISPOSAL AND RECYCLING OF THIS PRODUCT

Please note that this Adexi product is marked with this symbol:

This means that this product must not be disposed of together with ordinary household waste, as electrical and electronic waste must be disposed of separately.

In accordance with the WEEE directive, every member state must ensure correct collection, recovery, handling and recycling of electrical and electronic waste. Private households in the EU can take used equipment to special recycling stations free of charge. In some member states you can, in certain cases, return the used equipment to the retailer from whom you purchased it, if you are purchasing new equipment. Contact your retailer, distributor or the municipal authorities for further information on what you should do with electrical and electronic waste.

GUARANTEE TERMS

The guarantee does not apply:

- if the above instructions are not followed
- if the appliance has been interfered with
- if the appliance has been mishandled, subjected to rough treatment, or has suffered any other form of damage
- if faults have arisen as a result of faults in your electricity supply.


Due to the constant development of our products in terms of function and design, we reserve the right to make changes to the product without prior warning.

IMPORTER

Adexi Group

www.adexi.eu

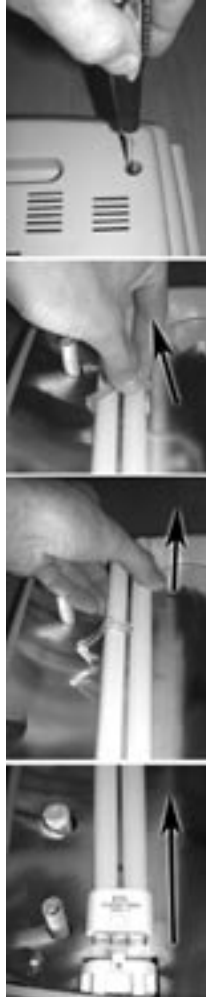We cannot be held responsible for any printing errors.

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