631-033

SE	Elektronisk personvåg med kropps fettanalys	2
DK	Elektronisk personvægt med kropsfedt-måler	7
UK	Electronic personal scale with body fat analysis.....	12
NO	Elektronisk personlig vekt med analyse av kroppsfett	16
FI	Elektroninen henkilövaaka, jossa rasvaprosentin mittaus....	21
DE	Elektronische Personenwaage mit Körperfettmesser	26
PL	Elektroniczna waga domowa z funkcj pomiaru procentowego udziaiu tkanki tłuszczowej	31
RU	Электронные весы для личного пользования с анализом содержания жира	35

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Der tages forbehold for trykfejl.

INTRODUCTION

Please read this instruction manual carefully and familiarise yourself with your new BIA scale before using for the first time. Please retain this manual for future reference.

SAFETY INSTRUCTIONS


- Read all instructions before using.
- Use this appliance only for its intended use as described in this manual.
- The readings of this body fat meter are only for reference. For the overweight people who want to reduce the body fat and underweight people who want to increase it, please consult your physician about medication and diet.
- People who are using pacemaker or other medical implant with electronic components, should not use the BIA scale since implants performance can be impaired
- In order to prevent interference: do not use the device near equipment with strong electrical field. Keep the device at a sufficient distance from such equipment or turn them off when using this scales.

SAVE THESE INSTRUCTIONS

SPECIAL SAFETY INSTRUCTIONS


- This device is suitable for adults and children above 5 years of age, not for people with pacemakers!
- The device will work correctly if your feet are bare and dry and if you are in standing position.
- The device has a maximum body weight limit of 150 kg. A larger weight than 150 kg may damage the device.
- The device may display incorrect results when measuring the body fat of pregnant women, adults over 70 years of age, dialysis patients or professional athletes.

INSERTING BATTERIES

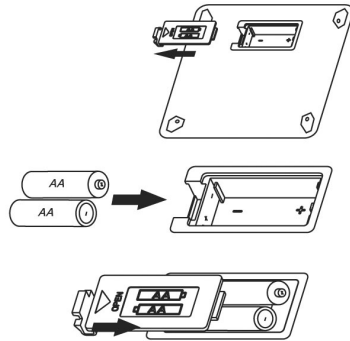Open the battery compartment on the back of the device.

Insert 2 pcs. AA type 1,5 V batteries in the battery holder as illustrated on the battery compartment cover.

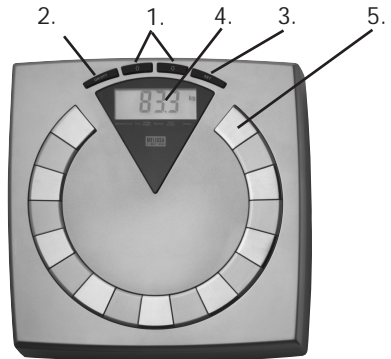
1. Arrows for selection
2. On/Off Button
3. "SET" button
4. Display
5. Contact plates


PROGRAMMING YOUR PERSONAL DATA

Attention: Carry out the following operations one after the other and within a few seconds or the device will automatically interrupt the programming mode.


1. Press the on switch to switch on the device "-HI-" appears briefly on the display followed by "1 - - 12".

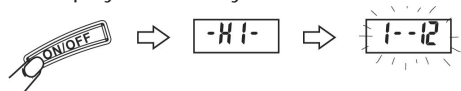
2. Press the "SET" key to enter programming mode.
3. Press the arrow keys to select the input unit mode (kg, cm / lbs, inch).
4. The "1" starts to flash on the display. Select a memory slot with the arrow keys and confirm your selection with the SET key. All the previous settings will be overwritten in the selected memory slot.

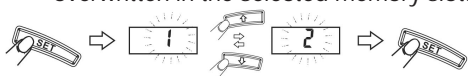5. Two symbols (male/female) appear on the display. Select your sex with the arrow keys and press the SET key.

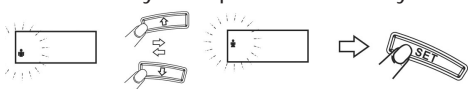6. Enter your height in cm. By pressing the arrow key again and press the SET key to confirm.

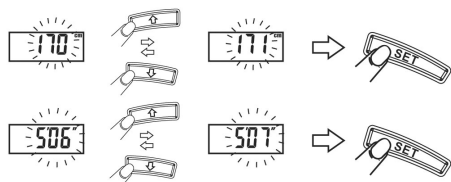7. The last value to be entered is your age. Use the arrow key again to select age, confirm by pressing the SET key.

8. "1 - - 12" appears in the display and shortly after "0,0". After a few seconds the device switches automatically off.

When age is 12 or below 12 the scale is in child mode.

WEIGHT AND FAT LEVEL DETERMINATION

First save your personal data in a memory slot.

1. Press the On switch to switch on the device. "1 - - 12" appears shortly in the display. Select the memory slot with your personal data
2. After a few seconds "0,0" appears in the display. Stand barefoot on the contact surface of the device
3. First your weight is shown on the display. Then the device determines your actual fat. The obtained value appears in percentage on the display (See fig. 1+2).
4. The device switches off automatically

Note: If "FFFF" appears in the display, your body weight exceeds 150 kg.

"FFFF". Weight overload
 "LO": Change battery
 "Err". Error

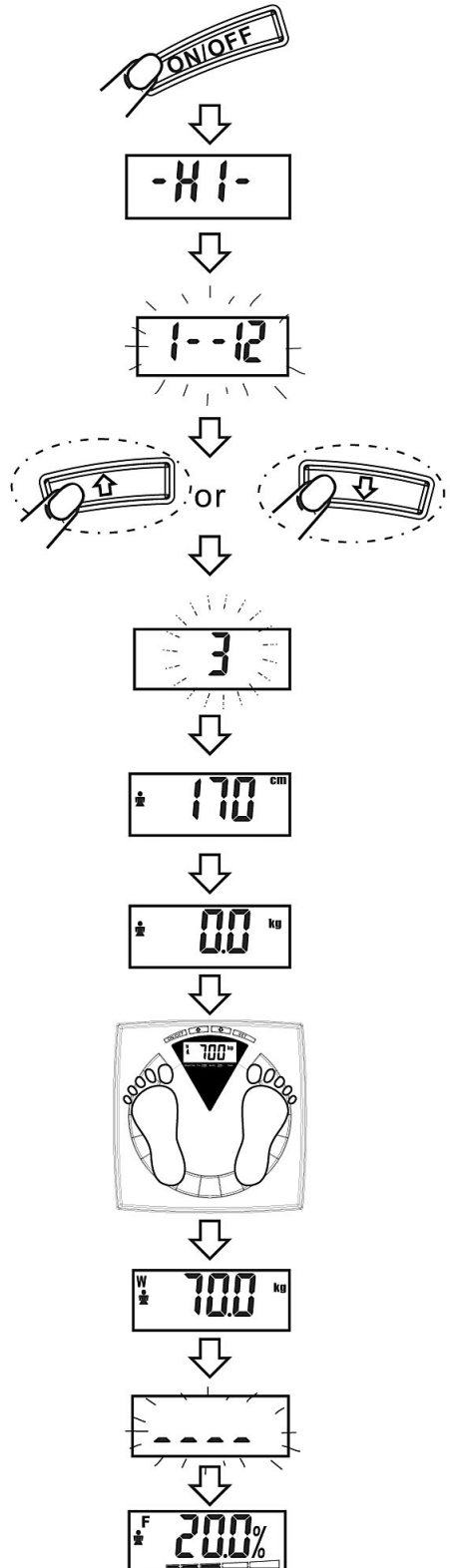
Fig. 1: Women


Age	Thin	Under weight	Normal	Over weight	Obese
6-12	3-7	8-13	14-21	22-29	30-38
13-17	3-8	9-14	15-23	24-30	31-38
18-39	3-13	14-20	21-25	26-31	32-38
40-55	3-14	15-21	22-26	27-32	33-39
56-75	3-15	16-23	24-28	29-34	35-41
Athlete			15-25		

Fig. 2: Men

Age	Thin	Under weight	Normal	Over weight	Obese
6-12	3-6	7-12	13-20	21-28	29-38
13-17	3-7	8-14	15-22	23-29	30-38
18-39	3-10	11-17	18-24	25-29	30-36
40-55	3-11	12-18	19-25	26-31	32-38
56-75	3-12	13-20	21-27	28-32	33-39
Athlete			3-15		

SIMPLE WEIGHING FUNCTION

1. Press the on switch (2) to switch on the device "HI" appears briefly on the display followed by "1 - - 12"
2. After a few seconds "0,0" will appear in the display. Stand on the contact surfaces of the device.
3. Your weight is shown on the display (3). After a few seconds the device shuts off automatically.

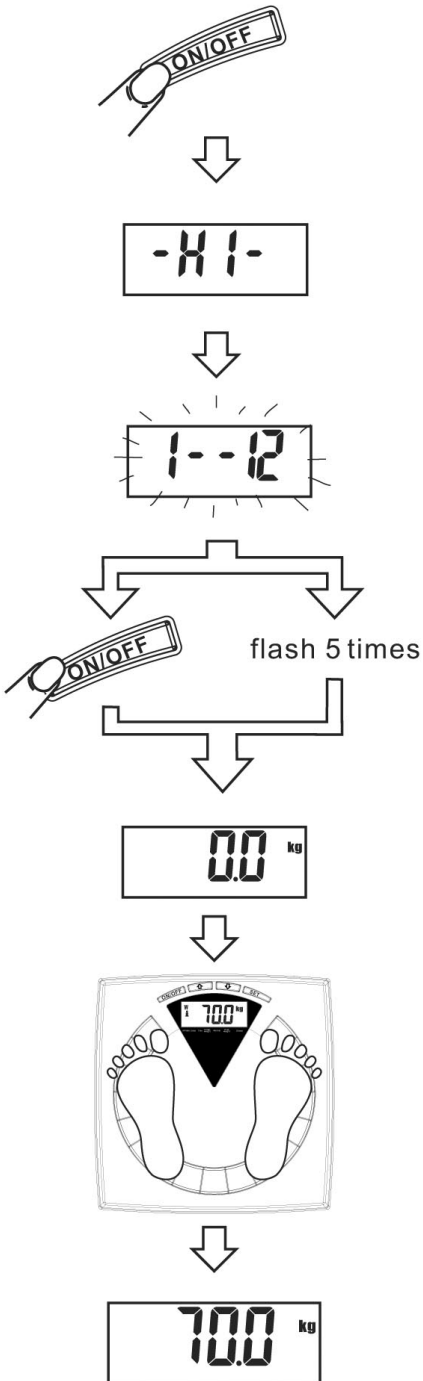MAINTENANCE:

1. Clean the appliance with a dry soft cloth. Do not use alcohol, benzene, thinner or other harsh chemicals to clean the device
2. Never immerse the appliance into water and other liquids!

ENVIRONMENTAL TIP:

An electronic appliance should, when it is no longer capable of functioning be disposed of with the least possible environmental damage according to the local regulations in your municipality. In most cases you can discharge the appliance at your local recycling centre.

THE WARRANTY DOES NOT COVER:

- If the above points have not been observed.
- If the appliance has not been properly maintained, if force has been used against it or if it has been damaged in any other way.
- Errors or faults owing to defects in the distribution system.
- If the appliance has been repaired or modified or changed in any way or by any person not properly authorised.

Owing to our constant development of our products on both functionality and design we reserve the right to change the product without preceding notice.

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Adexi A/S
Adexi AB