PHILIPS

sensiflex

dental center HX2745, HX1725

Download from www.Somanuals.com. All Manuals Search And Download.

ENGLISH	4
DEUTSCH	15
FRANÇAIS	27
NEDERLANDS	38
ESPAÑOL	49
ITALIANO	60
PORTUGUÊS	71


HX2745, 1725

Important

- ▶ Make sure the socket is live when you recharge the toothbrush. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is turned off.
- ▶ Charge the toothbrush for at least 16 hours before you use it for the first time.
- ▶ The cord of the appliance cannot be replaced. If the cord is damaged, you should have the appliance replaced by Philips or a service centre authorised by Philips in order to avoid hazardous situations.
- ▶ The mains plug of the appliance meets the latest European standards. The plug is safe, although the standard UK shaver socket holes may be slightly too wide for the pins of the plug. Make sure a proper connection is made.
- ▶ When you use the water jet for the first time or when you have not used it for a long time, it is necessary to set the volume control to its maximum position to make the water start flowing to the jet tip. If the water doesn't flow out after 10 seconds, put the jet tip in your mouth and suck it for a few seconds.
- ▶ If an "NiCd" indication is present on your appliance, the built-in rechargeable battery contains substances that may pollute the environment. When you discard the toothbrush in due course, remove the battery and dispose of it in an environmentally safe way.

Preparing the appliance for use

Mounting the base unit to the wall

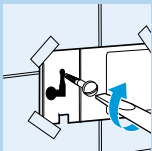
- 1** Mount the base unit to the wall, using 2 screws (3.5x30mm) and 2 wall plugs (S5) (not included).

To drill the holes in the right places, you can use the drilling stencil on the last page of these directions for use.

Do not allow the base unit to come into contact with water.

Using the toothbrush

Charging the toothbrush

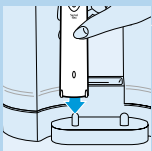- 1** Charge the appliance for at least 16 hours before using it for the first time.
Make sure the appliance is switched off when you start charging it.


► **Type HX2745 only:** when the battery is being charged, the green pilot light is on.

Tip: you can put the toothbrush in the base unit and let it charge continuously.


However, when you first start using the toothbrush, it is recommended that you use it until the battery is empty. To do so, use the toothbrush in the normal way without recharging it between brushing sessions. Recharge the

toothbrush when the motor almost stops running. Repeat this procedure 3 times in a row. To extend the lifetime of the battery, repeat this procedure every 6 months.


How to use the toothbrush

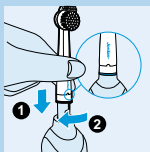- 1** Mount the brush attachment onto the toothbrush by inserting it into the bayonet catch and turning it clockwise a quarter of a turn.

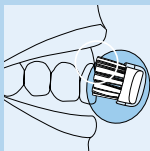- 2** Put toothpaste on the brush.

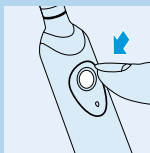- 3** Place the brush head in your mouth.

- 4** Switch the appliance on.

For the best brushing method, minimal brushing time and functions of the toothbrush, see the 'Brushing tips' section below.

Brushing tips


- 1** Place the brush head against the tooth in such a way that it covers the surface of the tooth as well as the gum line.
 - The round brush head removes plaque from the teeth and beneath the gum line.
 - The individually moving Active Tip cleans between the teeth.

2 Hold the brush head in this position for a few seconds and then move on to the next tooth.

- ▶ Make sure you clean all surfaces of the teeth (inside, outside and chewing surface) as well as the gum line properly. A fixed brushing pattern will prevent you from skipping certain areas.
- ▶ In general, 2 minutes will be enough to clean your teeth properly. Naturally, you can also brush longer than 2 minutes.

Type HX2745 only: a red pilot light will blink to indicate that 2 minutes have passed.


The Gum Protection system guarantees that you do not exert too much pressure on the teeth and gums.

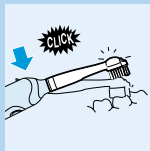
- ▶ If you exert too much pressure on the toothbrush, you will feel and hear a click when the brush head flexes backwards. As soon as you reduce the pressure, the brush head will return to its normal position. This will protect your gums.

Tip: if plaque is removed properly, the gums will have a pink colour and will firmly surround the teeth. Healthy gums do not bleed when you eat or brush your teeth. If, nevertheless, your gums bleed slightly the first few times you use the toothbrush, this bleeding will generally stop after a few days. If it does not, consult your dentist or dental hygienist.


Cleaning the toothbrush

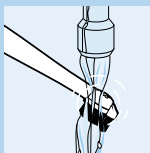
- 1 Switch the toothbrush on and rinse the brush attachment under the tap with lukewarm water (max. 60°C).

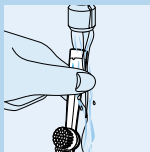- 2 Switch the toothbrush off. Detach the brush attachment and rinse the inside.

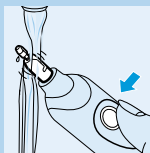- 3 Switch the toothbrush on again to rinse toothpaste residues from the drive shaft.

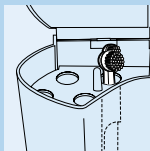- 4 Store the brush attachment in the brush compartment.


Brush head replacement

- 1 Replace the brush attachment every 3 months for the best brushing results.
 - Worn bristles will move plaque around rather than remove it. They may also damage the gums. It is therefore important that you replace the brush head every three months.

- ▶ *The brush attachment should also be replaced when the bristles start to bend outwards.*




You can obtain new brush attachments from the shop in which you purchased the appliance.


Using the water jet

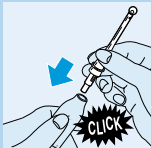When you use the water jet for the first time or when you have not been using it for a few days:

- 1** Click the jet tip onto the appliance.
- 2** Remove the water reservoir from the base unit.
- 3** Fill the water reservoir with lukewarm water.
- 4** Put the water reservoir back onto the base unit. Press it home properly.
- 5** Press the on/off button to switch the water jet on.

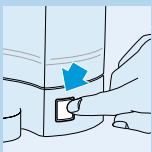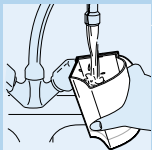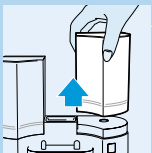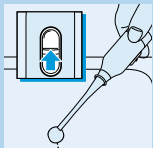
- 6** Push the slide on the handgrip upwards to make the water start running.

When you use the water jet for the first time, let the appliance run until the water reservoir is completely empty (let the water drain away into the washbasin). You can then refill the water reservoir and start using the appliance.


Tip: for extra freshness, you can add mouthwash to the water in the water reservoir. After using mouthwash, fill the water reservoir with clean water and let the appliance run until the water reservoir is empty (let the water drain away into the washbasin).

How to use the water jet

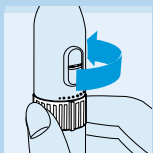
- 1** You can adjust the water volume by turning the control ring.


Tip: we advise you to use a low water volume setting while you are getting accustomed to the oral water jet and when you are troubled by gum inflammation.

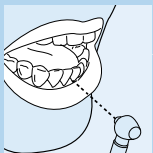
- 2** Bend over the washbasin and hold the jet tip close to your teeth.

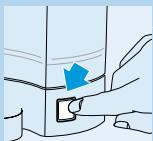- 3** Press the on/off button to switch the water jet on.

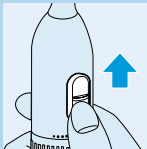- 4** Push the slide upwards to make the water start running.

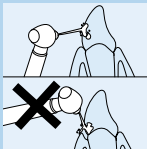
- 5** Guide the water jet along the front and the back sides of the teeth and along the gums. Direct the water jet at the spaces between the teeth in particular to remove any remaining food particles.


- ▶ *Make sure the water jet is at right angles to the teeth.*
- ▶ *Do not direct the water jet at the gums in such way that the jet penetrates beneath the gum line.*

- ▶ You can switch the water supply off with the slide before changing the position of the jet tip in your mouth.

After using the water jet

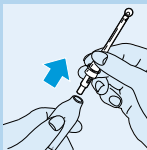- 1** Remove the jet tip and put it in the storage compartment.

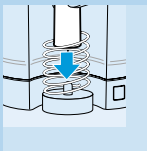- 2** Put the handgrip back into its holder.

Tip: you can obtain extra or replacement jet tips from your nearest Philips Service Centre.

Cleaning the water jet

- 1** Remove the mains plug from the wall socket.
- 2** Clean the appliance with a damp cloth and some detergent, if necessary.

Do not use abrasive or aggressive cleaning agents scourers, alcohol, acetone, etc. to clean the appliance.


- 3** Clean the storage compartment and the water reservoir in hot soapy water.
 - 4** Rinse them under a running tap.
- ▶ The storage compartment, the handgrip holder and the water reservoir may also be cleaned in a dishwasher.

Environment


When you discard the toothbrush in due course, remove the battery and dispose of it in an environmentally safe way.

- ▶ NB: you can also hand the appliance in at a Philips Service Centre, which will dispose of the battery for you.


The battery must be completely empty when you remove it. You should therefore let the appliance run until it stops before removing the battery.

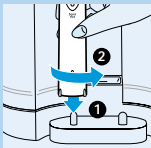
- 1** Place the toothbrush on the pin of the base unit. Turn the toothbrush anti-clockwise to unscrew the bottom of the toothbrush.

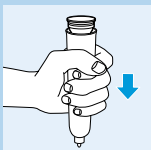- 2** Turn the toothbrush upside down and press it down on a hard surface to remove the battery holder.

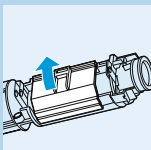- 3** Open the lid of the battery holder by means of a screwdriver.

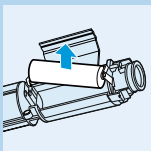- 4** Lift the battery out of the battery holder by means of a screwdriver.

After the bottom has been unscrewed, the appliance may not be watertight anymore and the guarantee will become invalid.

Information & service

If you need information or if you have a problem, please visit the Philips Web site at www.philips.com or contact the Philips Customer Care Centre in your country (you will find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, turn to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

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