TANITA

Monitoring Your Health

Digital Lithium


Operating Instructions – HD-332

OPERATING INSTRUCTIONS

Digital Lithium Scale – HD-332

- ❶ Measuring Platform
- ❷ LCD Display
- ❸ Battery Compartment on Bottom of Scale
- ❹ Weight Mode Switch on Bottom of Scale
- ❺ Arrow Keys

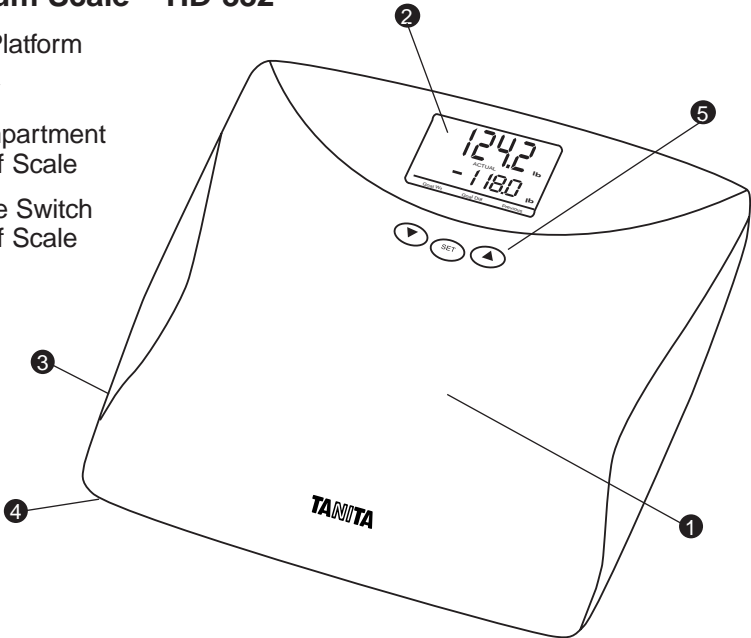Specifications

Maximum Capacity: 330 lb (150 kg)

Graduation: 0.2 lb (0.1 kg)

Power Supply: 1 x 3 volt CR2032 Replaceable Lithium Battery (Included)

UK Specifications

Maximum Capacity: 23 st 8 lb (150 kg) (330 lb)

Graduation: 1/4 lb (0.1 kg) (0.2 lb)

Power Supply: 1 x 3 volt CR2032 Replaceable Lithium Battery (Included)

Features

- Digital accuracy and reliability.
- Durable construction, non-skid surface.
- "L" message appears when battery needs replacing
- Auto Zero

Before Use

Select your preferred weight mode using the Weight Mode Switch

Placing Scale

The scale should be used on a hard even surface for accurate weighing. Not recommended for use on carpeted surfaces. Placing scale on a very soft, deep pile rug may affect weighing accuracy.

Goal Weight Input

This feature enables you to program your goal weight. Press “Set” to activate the scale for data input. The User Number “-1-” will flash on the first row of the LCD display. Press Up/Down arrow key to select desired User number (1 or 2) and press “Set” to confirm. The default Goal Weight will be displayed in the second row of the LCD display and there will be a triangle icon pointing to “Goal Wt.” Press the Up/ Down arrow keys to select your “Goal Weight” (press and hold the arrow key to advance quicker). Once your Goal Weight is reached, press “Set”. The target weight will flash to confirm that programming is complete. Scale will turn off automatically.

Goal Weight Function

After inputting your “Goal Weight”, this procedure will enable you to determine your “Goal Difference” by tracking changes between your “Goal Weight” and your “Actual Weight”. In addition, you will be able to recall your “Previous Weight”.

Start with scale off. Press any Arrow key until your User Number appears on the display. Press “Set” to confirm your number on the display. “00.0” will flash 3 times (Calibration Process). After 00.0 stops flashing, immediately step on scale and wait for weight reading to appear.

Note: “Err” will appear on the display if you stand on the scale while “00.0” is still flashing.

“Err” will appear on the display if you do not stand on the scale within 10 seconds after 00.0 stops flashing.

Readings

Top line will display “Actual Weight”. Bottom line will display “Goal Weight”, “Previous Weight” and “Goal Difference”. These 3 readings will flash every 3 seconds, a total of 2 times.

Note: “---” will flash at the bottom line if there is no target weight input from the selected User Number.

Weight Only Function

This procedure will enable you to determine your weight alone. You will not see “Goal Weight”, “Previous Weight” or “Goal Difference” in this weight only function.

Gently push down the middle of the platform with your foot and release the pressure to activate the scale. “CAL” will appear on the display. When “00.0” shows on the display, step on the scale (Fig. 1). Stand still while the display shows running 0000 for 3 seconds. The LCD will display your weight (Fig. 2). If you step off from the scale, the LCD will display your weight for 5 seconds then the scale will automatically shut off. If you remain on the scale, the LCD will display your weight for 10 seconds then automatically shut off.

FIG. 1

FIG. 2

Note: If you step on the scale without activating the scale, the display will show "Err" for 2 seconds then automatically shut off.

"Err" will appear on the display if you do not stand on the scale within 10 seconds after 00.0 stops flashing.

"Err" will appear on the display when the scale's weight capacity is exceeded.


Care and Maintenance

- Place the scale in an area free from heating equipment, high humidity, or extreme temperature changes.
- Never submerge the scale in water or spill chemical liquids on it. To clean, use a dampened cloth and mild detergent.
- Avoid excessive impact or vibration to the unit.

Replacing the Batteries

Lithium batteries can be typically used for many years before replacement is necessary.

When the low battery sign, L, appears on the LCD, please replace with one new lithium battery.

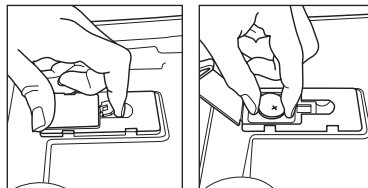Open the battery cover on the back of the measuring platform. Insert the new battery as indicated.

Note: Be sure that the polarity of the batteries is set properly. If you do not intend to use this unit for a long period of time, it is advisable to remove the batteries before storage.

After Sales Service

Please contact the Tanita office nearest you regarding operation or trouble, request for repair, and/or placement of an order for parts for products.

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Tanita Health Equipment H.K. Ltd.
Tanita Europe GmbH
Tanita UK Ltd.