## Precautions

1. Clean exposed parts with a soft, slightly damp cloth. To remove stains, use a mild soap. Never use detergents, excess water, treated cloths, harsh cleaning agents, or sprays. Do not attempt to lubricate or open scale casing as this will void your warranty.
2. Although your scale is designed to be maintenance-free, it contains sensitive electronic parts. Avoid rough treatment. Do not drop, kick, or jump on it. Treat it with care to ensure the best performance.
3. This scale is designed to read weights between 11 and 440 lbs / 5 and 200kgs / 2 and 31 st 6 lbs. To prevent damage, do not place more than 440lbs / 200 kgs /31 st 6 lbs on it.
4. Do not store scale where you store cleaning chemicals. The vapors from some household products may affect the electronic components of your scale. Do not store scale on its side.
5. This scale is an extremely sensitive weighing device. In order to prevent run down of the battery, do not store anything on the scale.
6. Do not dispose of batteries in fire. Batteries may explode or leak. Remove the battery if the scale will not be used for a long period of

## Five (5) Year Limited Warranty

This scale is warranted against defects in materials of workmanship five (5) years of the original purchaser from date of retail purchase. It does not cover damages or wear resulting from accident, misuse, abuse, commercial use, or unauthorized adjustment and/or repair. Do not return to retailer. Should this scale require service (or replacement at our option) while under warranty, please pack the item in the original packaging and return it prepaid, along with store receipt showing date of purchase and a note explaining reason for return to:

**Taylor Precision Products**  
**2220 Entrada Del Sol, Suite A**  
**Las Cruces, New Mexico 88001 USA**

There are no express warranties except as listed above. This warranty gives you specific legal rights, and you may have other rights which vary from state to state.

For additional product information, or warranty information outside the USA, please contact us through [www.taylorusa.com](http://www.taylorusa.com).

©2009 Taylor Precision Products and its affiliated companies, all rights reserved. Taylor® and Leading the Way in Accuracy® are a registered trademarks of Taylor Precision Products and its affiliated companies. All rights reserved.

**Made to our exact specifications in China.**  
**Not legal for trade.**  
**[www.taylorusa.com](http://www.taylorusa.com)**

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Leading the Way in Accuracy®

## Instruction Manual

Thank you for purchasing a Taylor® Lithium Electronic scale. All of our precision instruments for weight measurement are engineered and designed to meet the highest quality standards ...to assure you years of uncompromising accuracy and consistently dependable, convenient performance.

Please read this instruction manual carefully before use. Keep these instructions handy for future reference.

For more information on NBC's The Biggest Loser, or to join The Biggest Loser Club, please visit the official Biggest Loser website at [www.nbc.com/The\\_Biggest\\_Loser.com](http://www.nbc.com/The_Biggest_Loser.com) or [www.biggestloserclub.com](http://www.biggestloserclub.com).

## General Operation

1. There is a red strip under the Lithium battery cover to prevent battery drain in shipping. Please remove the strip before continuing with these instructions. Also, some models have a static cling label on the lens to prevent scratching. Please remove before use.
2. Your scale has been set at the factory to read weight in pounds. You may select kilograms or stone units of measurement by pressing the "kg/st/lb" button on the bottom of the scale.
3. Place scale on flat, hard surface. Carpeted or uneven floors may affect accuracy.
4. Step on the scale (disregarding the reading) and step off. This process insures the scale is set and is now ready for use. Do this whenever the scale is moved. **At all other times step straight onto the scale.**
5. Step onto the scale platform and remain still. The display will count up to your weight. Your weight is now displayed.
6. If you step off of the scale, the display will turn off automatically within 10 seconds. If you continue to stand on the scale, your weight will be displayed for 10 seconds and then the scale will turn off automatically.

## Battery Information

Your Taylor® scale has a lithium battery installed to provide extra long energy life for your scale. Under heavy usage (Health Club, Wellness Clinic, etc.) the battery may eventually need to be replaced.

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
## Battery Replacement

1. Use one 3 volt CR2032 lithium battery (as installed).
2. Open the battery compartment on the bottom of the scale by pushing the tab and lifting up the cover.
3. Remove the old battery from the compartment.
4. Place the new battery with the "+" sign up into the battery compartment.
5. Push gently down on the compartment cover until it clicks shut.
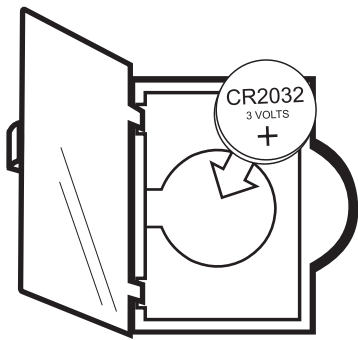
## Low Battery

Your Taylor® scale is equipped with a low battery indicator. It will display "LO" or readings will grow dim or become irregular if the battery needs replacement.

## Trouble Shooting

### In Case of Difficulty...

1. Place your feet in the same position every time you weigh yourself. Stand still during measurement.
2. Weigh yourself without clothing whenever possible. Clothes can vary in weight.
3. Move the scale to a location where it won't be bumped. If your scale is bumped while being activated or while in use, an inaccurate reading may occur.
4. For the most consistent readings, weigh yourself at the same time of day. Allow for weight fluctuations due to food or drink in your system.
5. The maximum capacity of this scale is 440 lbs / 200 kgs / 31 st 6 lbs. Should "----" appear on the display during a weighing, this indicates that the maximum weighing capacity of the scale has been exceeded. Remove the weight immediately; otherwise, permanent damage may occur.

## Reaching Target Weight

Maintaining proper weight is essential to healthy living. Obesity has been linked to a higher risk of heart disease, diabetes and some forms of cancer. The Body Mass Index (BMI) is the most widely accepted measurement of weight according to health professionals. Use the chart below to determine your own BMI by matching your height in the left hand column with your weight in the center. Anyone, male or female, with a BMI of 25 to 29 is considered over weight and a body mass index of 30 is considered obese. Remember to eat right and exercise.

Get a free diet profile at [www.biggestloserclub.com](http://www.biggestloserclub.com).

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