
# Waterpik® Water Flosser plus Sonic Toothbrush Model WP-900

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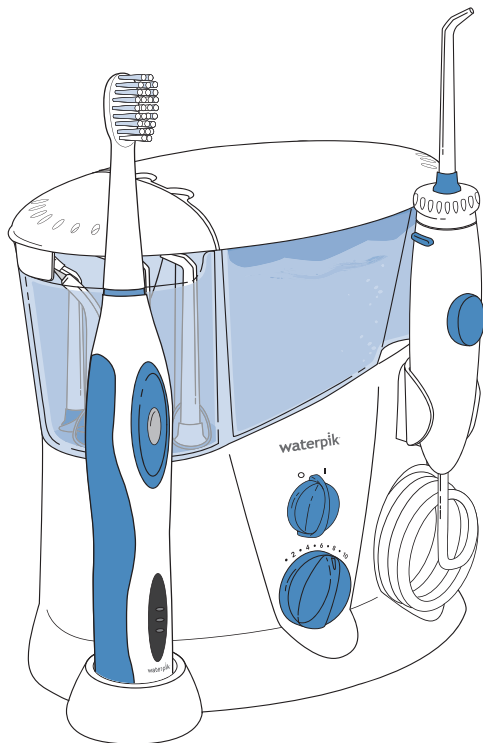.....  
[www.waterpik.com](http://www.waterpik.com)

Congratulations! By selecting the Waterpik® Complete Care you are now on your way to better gum health, and a healthier smile!

With a clinically proven combination of water pressure and pulsations, the Waterpik® Water Flosser removes harmful bacteria deep between teeth and below the gumline – where brushing alone won't reach.


We are so convinced that you will be happy with your new Waterpik® Complete Care, that we back it up with a 14-day guarantee. If you aren't satisfied that your gums are healthier after using the Waterpik® Complete Care for 14 consecutive days, feel free to return it to us for a full refund.

Thank you for choosing the Waterpik® Complete Care.

## Traditional Oral Care


Brushing alone can't easily reach hidden food debris and harmful bacteria.

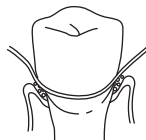Flossing has limited access below the gumline.

## Oral Care With The Waterpik® Water Flosser

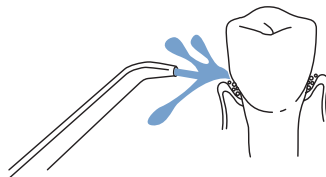The Waterpik® Water Flosser removes trapped food debris and harmful bacteria.

## TABLE OF CONTENTS

Important Safeguards	5
Product Description	7
Getting Started With Your Water Flosser	9
Care of Your Water Flosser	13
Water Flosser Troubleshooting Guide	16
Water Flosser Limited Two-Year Warranty	17
Getting Started With Your Sensonic® Professional Plus Toothbrush	18
Care of Your Sensonic® Professional Plus Toothbrush	20
Sensonic® Professional Plus Toothbrush Troubleshooting Guide	22
Sensonic® Professional Plus Toothbrush Limited Two-Year Warranty	23


## IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

### READ ALL INSTRUCTIONS BEFORE USING.

### DANGER: To reduce the risk of electrocution:

- Do not handle plug with wet hands.
- Do not immerse in water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into any liquid. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.
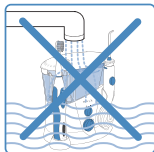
## WARNING:


To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. or your local distributor which can be found at [www.waterpik.com](http://www.waterpik.com).
- Do not direct water from the Water Flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.


- Do not use either the Water Flosser or Sonic Toothbrush where oxygen or aerosol sprays are being used.
- Do not use iodine, bleach, or water insoluble concentrated essential oils in the Water Flosser. Use of these can reduce product performance and will shorten the life of the product.
- Fill Water Flosser reservoir with water or other dental professional recommended solutions only.
- Do not use either the Water Flosser or Sonic Toothbrush while wearing any oral jewelry. Remove jewelry prior to use.
- Do not use either the Water Flosser or Sonic Toothbrush if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist and/or medical professional before using the Water Flosser or Sonic Toothbrush or any other oral hygiene aid.
- This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device containing magnets or electronics, contact your physician or the device manufacturer prior to using the Sonic Toothbrush.
- For household use only. Do not use the Water Flosser or Sonic Toothbrush outdoors.
- Closely instruct and supervise children and individuals with special needs in the proper use of the Water Flosser and Sonic Toothbrush.
- Do not remove the Pik Pocket™ flexible tip from the clear shaft. If it separates for any reason, discard both the tip and shaft and replace with a new Pik Pocket™ Tip. Pik Pocket™ Tip not included in all models.
- Do not use the Water Flosser for more than 5 minutes in each two-hour period. Longer usage can overheat the motor and cause product failure as well as void the warranty.
- Children should be supervised to ensure that they do not play with the appliance.
- This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug is intended to fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way.
- Be sure only the bristles of the toothbrush touch your teeth and gums when using the sonic toothbrush.
- This product does not contain any user serviceable parts in the toothbrush handle or Water Flosser unit. Do not attempt to open either component except for battery disposal. See page 21, "Battery Removal."
- Using toothpastes containing baking soda or peroxide, abrasive or chemical cleaners may damage the plastic housing of this product. Use of these products can shorten the life of this product.

**SAVE THESE INSTRUCTIONS.**


# PRODUCT DESCRIPTION

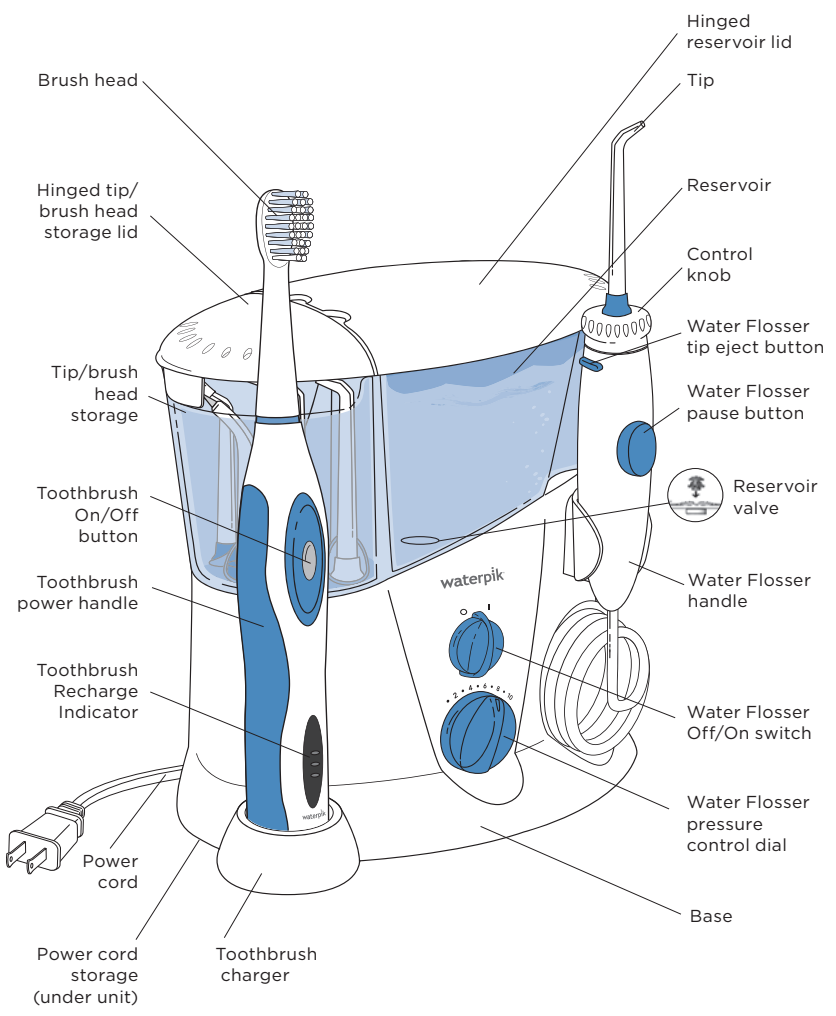
**NOTE:** all tips and brush heads not included in all models.

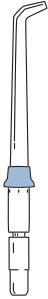
Classic Jet  
Tip

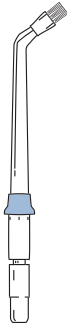Plaque Seeker®  
Tip

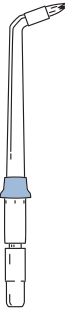Orthodontic  
Tip

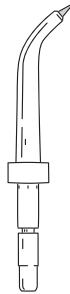Pik Pocket™  
Tip

Tongue  
Cleaner Tip

Sensonic®  
Standard  
Brush Head

Sensonic®  
Compact  
Brush Head

Replacement tips/attachments may be purchased directly from Water Pik, Inc. Call 1-800-525-2774 (USA) or 1-888-226-3042 (Canada) or visit our web site at [www.waterpik.com](http://www.waterpik.com).



## GETTING STARTED WITH YOUR WATER FLOSSER

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.

### Preparing the Reservoir

Lift the reservoir lid, remove the reservoir, and fill with lukewarm water. Replace the reservoir over the base and **press down firmly** and close the lid.

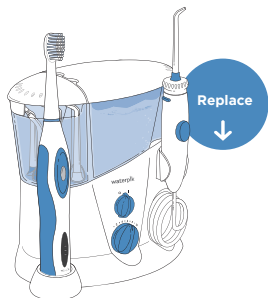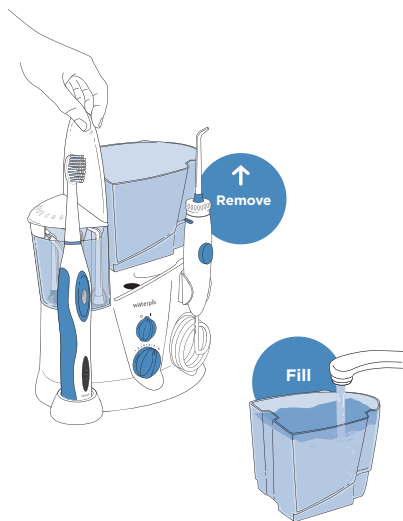
### Inserting and Removing Tips

Insert tip into the center of the Water Flosser handle, pressing down firmly until tip clicks into place. The colored ring will be flush with the top of the handle if the tip is correctly locked into place. Gently pull on tip to confirm it is seated properly and locked in place.


To remove tip from the handle push the tip eject button on the handle. The control knob will rise, releasing the tip. Remove the tip.

**Caution:** Do not press tip eject button while unit is running.

### Tip Insertion

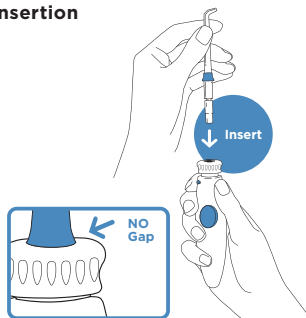
### Tip Removal

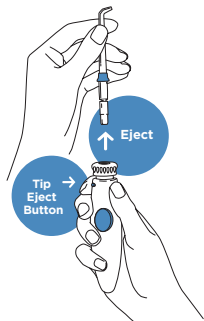
### Prime Unit (First Time Use Only)

Prime unit prior to very first use of the unit: Set pressure control dial on base to high (pressure setting No. 10), point the handle and tip into sink. Turn unit ON until water flows from the tip. Turn unit OFF.

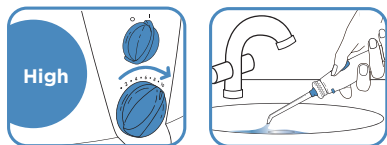
### Adjusting the Pressure Setting

Turn the pressure control dial on the base of the unit to the lowest setting (pressure setting No. 1) for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

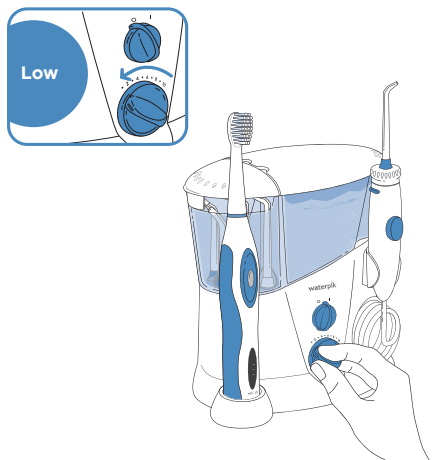
### Recommended Technique

Direct the stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue process until you have cleaned the inside and outside of both the upper and lower teeth.

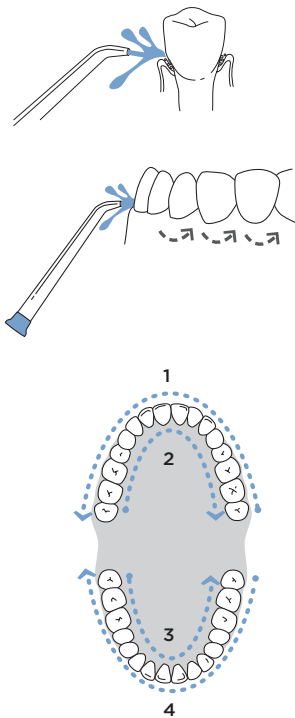
## Directing Water Stream

Direct the stream by turning the knob at the top of the handle.
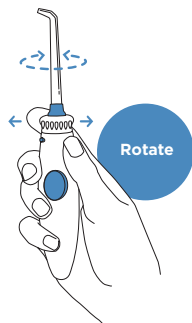
## Pause Control

You can temporarily stop the flow any time by pressing the pause button on the handle.

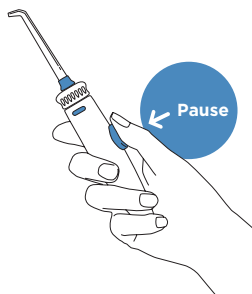
## Placement in Mouth

Lean low over sink and place the tip in mouth. Aim the tip toward teeth. Set pressure at desired setting. With your free hand, turn the unit ON.

## When Finished

Turn the unit OFF. Empty any liquid left in the reservoir.

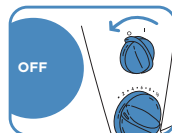
**Note:** If you used a mouthwash or other solution besides water, see "Using Mouthwash or Other Solutions" section.

## TIP USAGE

### Pik Pocket™ Tip

The Pik Pocket™ Tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets. It is not intended for general use.

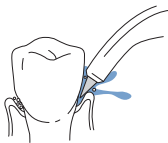To use the Pik Pocket™ Tip, **set the unit to the lowest pressure**

**setting.** Lean low over your sink and place the soft tip against a tooth at a 45-degree angle


facing the gum line and gently move the tip until it is slightly below the gum, into the pocket. Turn the unit on and continue tracing along the gum line, pausing briefly between the teeth until all periodontal pockets are treated.

**NOTE:** This tip must be used only with the Pressure Control set to 1. Use on a higher setting will shorten life of tip and Water Flosser handle.


### Tongue Cleaner Tip

To use the tongue cleaner tip, set the unit to the lowest pressure setting. Gently place the tongue cleaner in the center/midline of your tongue, about halfway back and start the Water Flosser. With light pressure pull forward. You can increase pressure as you prefer.
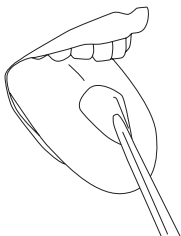
### Orthodontic Tip

The special Orthodontic Tip is uniquely designed to simultaneously brush and rinse hard-to-reach areas around braces and other dental work. Start on the last molar (back tooth). Aim water stream at 90 degrees to the tooth at the gumline. Gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth. Continue process until you have cleaned the inside and outside of both the upper and lower teeth.
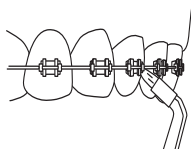
### Plaque Seeker® Tip


Direct the jet stream at a 90-degree angle to your gumline. Close lips slightly to avoid splashing but allow water to flow freely from mouth into sink. Start in the molar (back teeth) area working toward the front teeth. Place the tip close to the teeth so that the bristles are gently touching the teeth.

Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth. Continue until all areas around and between the teeth both inside and outside, top and bottom have been cleaned.


## Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. Ask your dental professional to recommend the appropriate solution for you. **NOTE:** use of some antibacterial solutions could shorten the life of your Waterpik® Water Flosser, see warnings section for details.

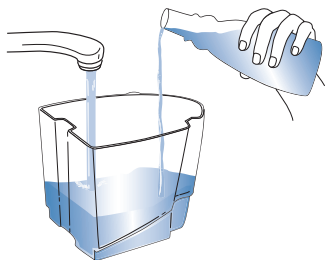After using any special solution in the reservoir, take a moment after every use to rinse the unit. Partially fill the reservoir with warm tap water. Aim the jet tip into your sink. Turn on the unit and let it run until empty. This will help keep the tip from clogging and help your unit to last longer.
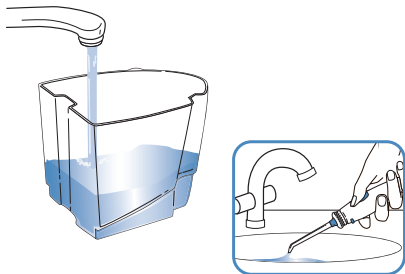
## Cleaning Your Waterpik® Water Flosser

Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

The tip storage and the water reservoir are top rack dishwasher safe.


**NOTE:** Remove reservoir valve before placing reservoir in dishwasher.

Remove


Replace

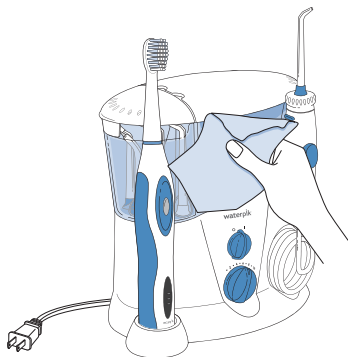Before exposing the Water Flosser to freezing temperatures, remove the reservoir and run the unit until completely empty.

## Accessory Storage


The accessory storage container has a hinged lid, which allows for easy access to your tips and brush heads. The storage container will hold both Water Flosser tips and Sensonic® brush heads. The storage container is also removable for easy cleaning.

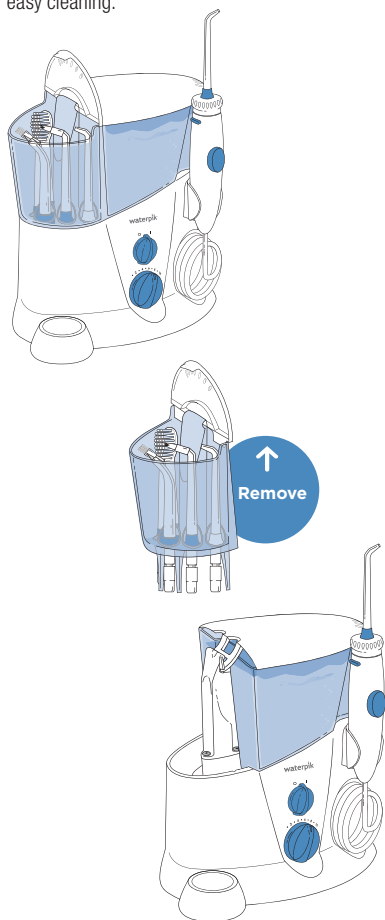
## Excess Power Cord Storage

The Water Flosser includes an area for excess power cord storage. After determining the length of power cord necessary you may push the extra cord into the cavity on the bottom of the unit.


First remove the cord from the retention slots.

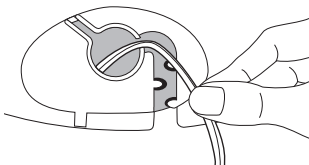Push the excess cord into the hole on the bottom of the unit until the desired length of cord is left.

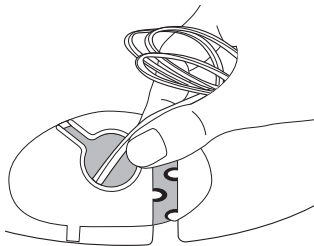Re-insert the cord into the retention slots.

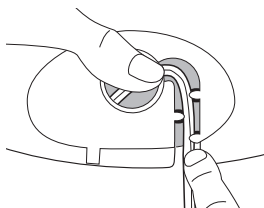
## Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all accessories go to [www.waterpik.com](http://www.waterpik.com).

You may also call our toll-free customer service number at 1-800-525-2774 (USA) or 1-888-226-3042 (Canada).

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

## RECOMMENDED TIP AND BRUSH HEAD REPLACEMENT INTERVALS

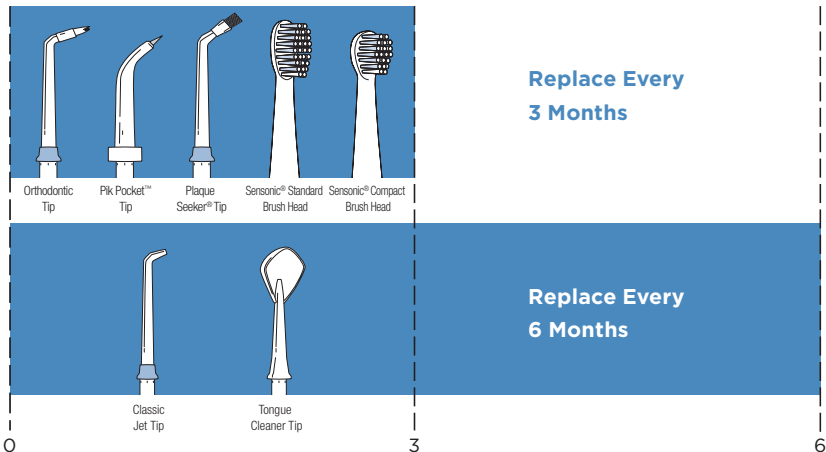Replacement tips/attachments may be purchased directly from Water Pik, Inc. Call 1-800-525-2774 (USA) or 1-888-226-3042 (Canada) or visit our web site at [www.waterpik.com](http://www.waterpik.com).

Download from [Www.Somanuals.com](http://Www.Somanuals.com). All Manuals Search And Download.

## WATER FLOSSER TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	CORRECTIVE ACTION
<b>Tips</b> Leakage between the tip and handle knob.	1) Tip not fully engaged.	1) Remove tip and reinstall.
<b>Pik Pocket™ Tip</b> Tear in soft tip.	1) Tip is worn. 2) Tip is split.	1) Replace with new tip. 2) Replace with new tip and use ONLY on lowest pressure.
<b>Water Flosser Base</b> Hose to the handle leaks.	1) Damage to the hose.	1) Visit <a href="http://www.waterpik.com">www.waterpik.com</a> or contact Water Pik, Inc. customer service for a self service repair kit.
Reservoir leaks.	1) Missing black reservoir valve.  2) Black reservoir valve is upside down.	1) Replace reservoir valve. Contact Water Pik, Inc. customer service for a replacement valve. 2) Ensure valve is right side up (see image to right). 
Inadequate pressure.	1) Reservoir not positioned properly.  2) Tip is clogged with debris. 3) Black reservoir valve is clogged.  4) Hard water deposits.	1) Re-position reservoir. Press down firmly on the reservoir. 2) Back flush tip under faucet or replace tip. 3) Remove reservoir valve and clean by massaging the valve under warm running water and place back into reservoir. Ensure the valve is right-side up (see above image). 4) Add 2 tablespoons of white vinegar to a full reservoir of warm water. Set pressure control dial to high, point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.
Unit doesn't start.	1) Motor not running.	1) Be sure outlet is functioning; test with another appliance (check electrical reset button and/or wall switch).

Still have questions? Go to [www.waterpik.com](http://www.waterpik.com) or call Water Pik, Inc. Customer Service.  
 Toll -Free Number: 1-800-525-2774 (U.S.A.) 1-888-226-3042 (Canada).

Download from [Www.Somanuals.com](http://Www.Somanuals.com). All Manuals Search And Download.

## WATER FLOSSER LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips and brush heads.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

## GETTING STARTED WITH YOUR SENSONIC® PROFESSIONAL PLUS TOOTHBRUSH

### Charging the Sensonic® Professional Plus Toothbrush


#### Charge for 24 hours before first use.

The charge indicator lights on the handle will go through a fill-up sequence when charging and all 3 lights will stay on continuously when the unit is fully charged.

A full charge will last for approximately 1 week of typical use (2 minutes of brushing twice a day).



24 Hours

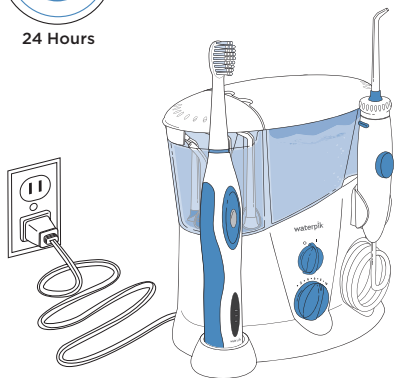If the toothbrush is to be stored for an extended period of time (six months or more), charge it for 24 hours before storing. Never store the toothbrush for extended periods unless it is fully charged.

### Recharge Indicator Lights

When the unit is taken off a powered charging base it will indicate the level of battery charge for approximately 10 seconds and then the lights will go off. The lights will stay illuminated for 10 seconds after use indicating the level of battery charge.

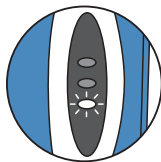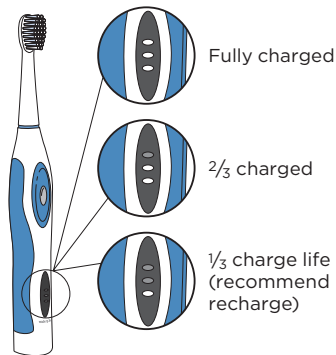
Unit will not operate if the bottom LED begins to blink rapidly.


In order for the toothbrush to begin to work again, it must be placed on the charger for approximately 1 hour. It will still require a 24 hour charge to regain full charge status.

## Sonic brushing



- Slide the brush head onto handle and push it down firmly.
- Wet the bristles and apply a small amount of toothpaste.
- Place the brush head into your mouth before turning the unit on.
- Select low or high brushing speed by pressing power button. Press button once for low speed, and twice for high speed.
- If you hold the on/off button down for approximately 1.5 seconds while in the low speed position the brush will shut off.
  - Low speed – gentle cleaning and gum massage
  - High speed – superior plaque removal
  - For best results, position brush along gumline at a 45-degree angle.

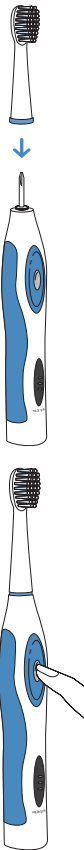- Use light pressure with a small back and forth motion covering one or two teeth at a time. Overlap the teeth as you work your way around the mouth. Brush first the outsides, then the insides, finally the chewing surfaces of the teeth.
- Pause briefly between the teeth to allow the longer bristles to clean between your teeth.
- Move the brush slowly around entire mouth. Brush each quadrant of your mouth thoroughly for 30 seconds.
- **Quad Timer:** Your Sensonic® Professional Plus Toothbrush will create a short pause every 30 seconds, to remind you to move to the next quadrant of your mouth.

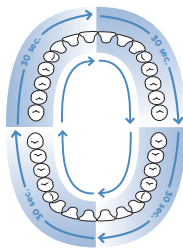- **2-Minute timer:** The brush will automatically shut off after 2 minutes of run time.

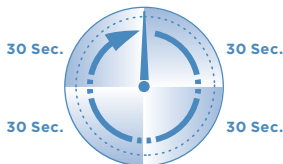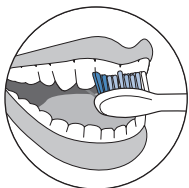
**2 Minutes**

## CARE OF YOUR SENSONIC® PROFESSIONAL PLUS TOOTHBRUSH

### Using the interdental brush


(Not included in all models)

- Slide the brush head onto handle and push it down firmly.
- Place the brush head into your mouth before turning the unit on.
- Brush may be used with or without toothpaste.
- Place brush tip at a 90-degree angle to the gumline between two teeth.
- Turn on power handle and lightly follow crevices between teeth.

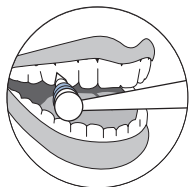
### Cleaning

1. After each use, turn off unit, remove the brush head from the handle and rinse brush head under running water.

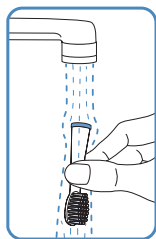2. Wipe the handle dry after each use.
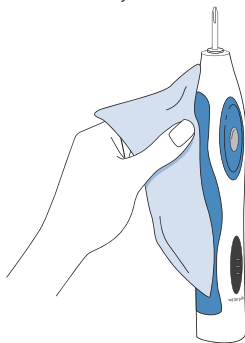
### Service

Waterpik® toothbrushes have no consumer servicable or replaceable electrical items and do not require routine service maintenance. For all accessories go to [www.waterpik.com](http://www.waterpik.com).

You may also call our toll-free customer service number at 1-800-525-2774 (USA) or 1-888-226-3042 (Canada).

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Replacement tips/attachments may be purchased directly from Water Pik, Inc. Call 1-800-525-2774 (USA) or 1-888-226-3042 (Canada) or visit our web site at [www.waterpik.com](http://www.waterpik.com).

### Disposal

This product contains rechargeable batteries. The batteries are not replaceable.

Please do not dispose of the unit in the household waste. Remove the batteries per the instructions on page 21 and recycle at your local recycling center. In so doing, you will help preserve the environment.



## Battery Removal

### Note: This process is not reversible



The rechargeable battery inside the toothbrush cannot be replaced but can be removed for recycling. Remove the battery only at the end of the product's life. Doing so earlier will destroy the product. Deplete the battery charge by running 2 minute cycles until the brush no longer operates. Observe basic safety precautions.

At the end of the product's useful life, open the handle as shown, remove the batteries and dispose of them according to local regulations. The rest of the product can be discarded per local regulations.

1. Use coin or washer to rotate cap counter-clockwise until cap protrudes. Rotation will stop at this point.
3. Remove cap and pull internal body out. Cut tape and remove.

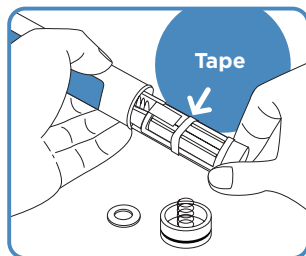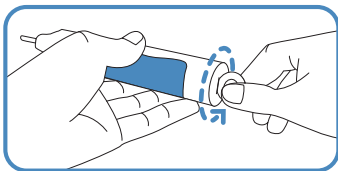2. Push down hard with shaft on hard surface.
4. Lift battery pack up and cut battery wires, one at a time, to free from circuit board.
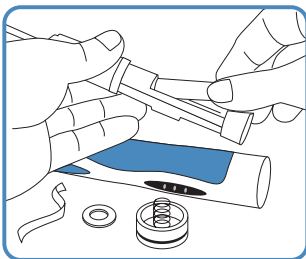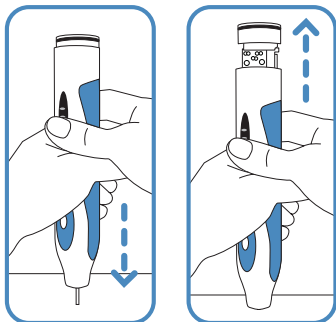
## SENSONIC® PROFESSIONAL PLUS TOOTHBRUSH TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	CORRECTIVE ACTION
I placed a discharged toothbrush on the charger, but nothing happens.	1) Plug may not be properly inserted into the outlet. 2) Confirm outlet is live.	1) Insert plug properly into outlet. 2) Test the outlet with another appliance.
I placed a completely discharged toothbrush on the charger, but only one light flashes and unit still does not turn on.	1) A fully discharged power handle will not operate until it has been charged for at least 1 hour.	1) Charge for at least 1 hour. The power handle will require a 24 hour charge to be fully charged.
I had not used this product for several months and when I tried to use it the handle would not turn on. <b>OR</b> After several months of no use, the battery is dead.	1) A battery-operated device not used for a long time will gradually lose charge. The power handle will gradually lose its charge.	1) Charge the handle for 24 hours. Charge the power handle fully prior to storage.
While using this toothbrush, I see the lights flash and the vibration stops momentarily.	1) The toothbrush has a built-in timer allowing vibration to pause and the lights to flash every 30 seconds. The quadrant timer guides the user to move to the next quadrant of the mouth.	1) None. This is normal operation.
While I am using the toothbrush it shuts off by itself.	1) The toothbrush has a timer that shuts off after 2 minutes of use, the recommended brushing time.	1) None. This is normal operation.

Still have questions? Go to [www.waterpik.com](http://www.waterpik.com) or call Water Pik, Inc. Customer Service.  
Toll-Free Number: 1-800-525-2774 (U.S.A.) 1-888-226-3042 (Canada).

Download from [Www.Somanuals.com](http://Www.Somanuals.com). All Manuals Search And Download.

Water Pik, Inc. warrants to the original consumer of this product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as brush heads.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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# waterpik®

## Water Pik, Inc.

1730 East Prospect Road  
Fort Collins, CO 80553-0001 USA

[www.waterpik.com](http://www.waterpik.com)

Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Benelux, Canada, Chile, China, Hong Kong, India, Japan, Korea, Mexico, New Zealand, Norway, Pakistan, Russian Federation, Singapore, South Africa, Switzerland, Taiwan, Ukraine, and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Australia, Brazil, Canada, China, EU, Hong Kong, Japan, Mexico, Korea, Russian Federation, Switzerland, Trinidad & Tobago, and the United States.

Sensonic® is a trademark of Water Pik, Inc. registered in Canada, China, EU, Japan, Switzerland, and the United States.

Treat Yourself Better™ and Pik Pocket™ are trademarks of Water Pik, Inc.

Plaque Seeker® is a trademark of Water Pik, Inc. registered in Canada, China, EU, and the United States.

**waterflosser®** is a trademark of Water Pik, Inc. registered in Canada, China, Mexico, and Switzerland.

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©2012 Water Pik, Inc.
Treat yourself better™ et Pik Pocket™ sont des marques déposées de Water Pik, Inc.
Treat Yourself Better™ y Pik Pocket™ son marcas de comercio de Water Pik, Inc.
**waterflosser®** es una marca de comercio de Water Pik, Inc., registrada en Canadá, China, México y Suiza.