# TOPCOM®

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## 1 Intended Purpose

The MM-1000 is a handheld massager that uses vibration as a means of alleviating muscular strains and pains, as well as promoting relaxation. It can be used to stimulate meridians and enhance blood circulation. A wonderful way to relieve stress.

## 2 Safety advice

**The device is only intended for the purpose described in this user guide. The manufacturer cannot be held liable for personal injuries or damages occurring as a result of inappropriate or careless use.**

### 2.1 General

- Always read the safety instructions carefully and keep this user guide for future reference.
- Follow the basic safety precautions for all electronic equipment when using this product.
- This product is not a toy. Don't let children play with it.
- The device is only for private use and is not intended for medical or commercial purposes. If you have any concern regarding your health, consult your doctor before using this product.
- This unit is not waterproof. Don't expose this product to rain or moisture.
- Make sure that the appliance is not too close to heat sources such as ovens, radiators or direct sunlight.
- Do not use the product if it is not working properly, dropped or damaged.
- Repairs to electrical appliances must only be carried out by qualified persons. Incorrect repairs can lead to considerable danger for the user. In the event of repairs, please contact our customer service or an authorized dealer.
- Don't disassemble this product. It doesn't contain parts that can be repaired by the user.

- Do not drop or insert object into any opening.
- Do not place the basic unit in a damp room and at a distance of less than 1.5 m away from a water source.
- Do not use the product for more than 15 minutes at a time. Extensive use could lead to the product's excessive heating and shorter lifetime.
- Don't disassemble this product. It doesn't contain parts that can be repaired by the user.

## **2.2 Batteries**

- Batteries can cause fire or burns when the terminals are short-circuit.
- Never throw batteries in fire as they may explode.
- Keep batteries away from small children.
- Leaking acid from batteries may damage the device permanently. Remove the batteries if the device is not going to be used for a long period.
- Handle damaged or leaking batteries very carefully.


## **2.3 Massage precautions**

- It is not allowed to use the massager:
  - On young children under the age of 16.
  - On persons suffering from a herniated disk or from abnormal changes to the vertebra.
  - Not in case of cardiovascular diseases and headaches.
  - Not in the case of swelling, burns, inflammation, open wounds, fresh wounds, hematomas, bleeding or other injuries to the back or legs.
  - During pregnancy.
  - In bed.
  - After taking medication or consuming alcohol (reduced levels of awareness).
  - On a person with insensitive skin or on a person with poor blood circulation.
  - In a car.
  - While operating machinery.

- If you are fitted with a cardiac pacemaker, implant or other aid.
- If you suffer from any of the following: common cold accompanied by fever, varicose veins, thrombosis, phlebitis, jaundice, diabetes, nervous system diseases (e.g. sciatica), tuberculosis, tumours, haemorrhoids, bruises or acute inflammation.
- Not in case of pains of unclarified origin.
- On the head.
- Close supervision is necessary when the appliance is used by or near invalids or disabled persons.
- Individuals with pacemaker should consult their doctor before use.
- Do not use this product before sleeping. The massage has a stimulating effect and can delay sleep.
- When the unit is working, do not put hands between the massage head and body.
- Only gentle force should be exerted against the unit in order to eliminate any/all risk of injury.
- If the material covering the massage mechanism becomes torn or damaged, discontinue use.
- Never put any body part (fingers or hands) inside the interior of the unit.
- Consult your physician, if you have any therapeutic questions regarding use of massager.
- Massage should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.

### 3 Button description

1. On/Off button
2. Battery cover
3. Illuminated massage head
4. Battery cover bolt

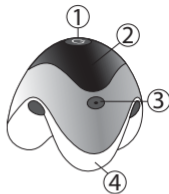
### 4 Getting started.

#### 4.1 Important

- Consult a licensed health care professional. He or she will advise you on whether it is safe for you to use a massager, what to do, and whether there are other things that you should be doing as well to help.
- Do not use directly on your skin. Use clothing or a cloth to protect your skin and aid in hygiene.
- Do not massage a body part more often than every second day.

#### 4.2 Installing batteries

- Press the battery cover bolt to release the battery cover. Watch very carefully how the battery cover is removed.
- Remove the battery security clip by pressing both ends.
- Install 3 AAA (1,5) batteries following the polarity as indicated in the battery compartment and re-install carefully the battery security clip and battery cover.
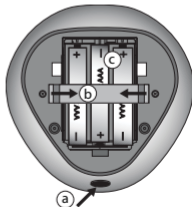
### 4.3 Using the massager


After installing batteries:

- Press the On/Off button (1) to activate the vibration. The massage heads will illuminate.
- Press the On/Off button (1) again to deactivate the vibration.

### 4.4 How and where to massage

- Keep the massager still on one of the next spots for approximately 30 seconds then move it onto the next part that needs massaging.

1. Neck
2. Back
3. Waist
4. Palm
5. Sole
6. Calf
7. Thigh
8. Arm
9. Neck and shoulder

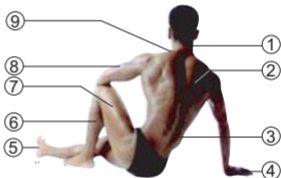- Apply little pressure at first.
- It is very important that the part you are massaging is relaxed. When you hold a massager, the shoulder on the side you hold the massager with will tighten. Therefore you should hold the massager with the opposite hand.
- Never press too hard.

### 4.5 Technical data

Dimensions (W/H/D)	90 x 90 x 70 mm
Weight	147 g (without batteries)
Power supply	3 x AAA (1,5V)

## 5 Cleaning

- Clean it with a soft cloth or slightly damp sponge.
- Never allow that water or any other liquid enters the device or accessories.
- Never use abrasive cleaners, brushes, gasoline, kerosene, glass polish or paint thinner to clean.

## 6 Disposal of the device (environment)

At the end of the product lifecycle, you should not throw this product into the normal household garbage but bring the product to a collection point for the recycling of electrical and electronic equipments. The symbol on the product, user guide and/or box indicate this.

Some of the product materials can be re-used if you bring them to a recycling point. By re-using some parts or raw materials from used products you make an important contribution to the protection of the environment. Please contact your local authorities in case you need more information on the collection points in your area.

## 7 Topcom warranty

### 7.1 Warranty period

The Topcom units have a 24-month warranty period. The warranty period starts on the day the new unit is purchased.

Consumables or defects causing a negligible effect on operation or value of the equipment are not covered.

The warranty has to be proven by presentation of the original purchase receipt, on which the date of purchase and the unit-model are indicated.

## 7.2 Warranty handling

A faulty unit needs to be returned to a Topcom service centre including a valid purchase note. If the unit develops a fault during the warranty period, Topcom or its officially appointed service centre will repair any defects caused by material or manufacturing faults free of charge.

Topcom will at its discretion fulfill its warranty obligations by either repairing or exchanging the faulty units or parts of the faulty units. In case of replacement, colour and model can be different from the original purchased unit.

The initial purchase date shall determine the start of the warranty period. The warranty period is not extended if the unit is exchanged or repaired by Topcom or its appointed service centres.

## 7.3 Warranty exclusions

Damage or defects caused by incorrect treatment or operation and damage resulting from use of non-original parts or accessories not recommended by Topcom are not covered by the warranty.

The warranty does not cover damage caused by outside factors, such as lightning, water and fire, nor any damage caused during transportation.

No warranty can be claimed if the serial number on the units has been changed, removed or rendered illegible.

**This product is in compliance with the essential requirements and other relevant provisions of the directive 2004/108/EC.**

**The Declaration of conformity can be found on:**

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