TOPCOM[®]

SMC-2000H Massage cushion

**USER GUIDE / HANDLEIDING / MANUEL D'UTILISATEUR
BEDIENUNGSANLEITUNG / MANUAL DE USUARIO
BRUKSANVISNING / BRUGERVEJLEDNING
ANVÄNDARHANDBOK / KÄYTTÖOHJE / MANUALE D'USO
MANUAL DO UTILIZADOR / UŽIVATELSKÁ PŘÍRUČKA
ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ / HASZNÁLATI ÚTMUTATÓ / INSTRUKCJA OBSŁUGI
MANUAL DE INSTALARE SI UTILIZARE / UŽIVATEĽSKÝ MANUÁL**

V.1.0 - 06/08

- UK** *The features described in this manual are published with reservation to modifications.*
- NL** *De in deze handleiding beschreven mogelijkheden worden gepubliceerd onder voorbehoud van wijzigingen.*
- FR** *Les possibilités décrites dans ce manuel sont publiées sous réserve de modifications.*
- DE** *Die in dieser Bedienungsanleitung beschriebenen Leistungsmerkmale werden unter dem Vorbehalt von Änderungen veröffentlicht.*
- ES** *Las características descritas en este manual pueden ser objeto de futuras modificaciones.*
- SE** *Funktionerna i denna bruksanvisning publiceras med reservation för ändringar.*
- DK** *Vi forbeholder os retten til ændringer af de specifikationer, der er beskrevet i denne brugsanvisning.*
- NO** *Funksjoner beskrevet i denne manualen kan endres uten nærmere informasjon.*
- FI** *Tässä ohjekirjassa julkaissut tiedot voivat muuttua ilman ennakoilmoitusta.*
- IT** *Le caratteristiche descritte nel presente manuale vengono pubblicate con riserva di modifica.*
- PT** *As características descritas neste manual são publicadas sob reserva de modificação.*
- CZ** *Možnost úpravy funkcí popsaných v této příručce vyhrazena.*
- GR** *Οι λειτουργίες που περιγράφονται στο παρόν εγχειρίδιο, δημοσιεύονται με επιφύλαξη τυχόν τροποποιήσεων.*
- HG** *A jelen használati útmutatóban leírt adatokra a gyártó fenntartja a változtatás jogát.*
- PL** *Właściwości opisane w niniejszej instrukcji obsługi są publikowane z zastrzeżeniem prawa wprowadzenia zmian.*
- RO** *Descrierile din acest manual sunt publicate sub rezerva modificarilor.*
- SK** *Vlastnosti popísané v tejto príručke sú publikované s vyhradeným právom na zmenu.*

1 Intended Purpose

The SMC-2000H is a portable shiatsu massage cushion that contains internal electronic motors and gears designed to massage the person sitting on it.

The full intensity Shiatsu function imitates a strong, stimulating massage.

The massage cushion can be fixed on an office chair or any other suitable seating with backrest (chair, sofa, or similar).

With the controller you can vary the type, location, or intensity of massage.

The infrared heating improves the comfort during winter.

2 What is Shiatsu

The translation of the Japanese word 'Shiatsu' is 'finger pressure'.

The primary aim of Shiatsu is to promote energy flow, which serves to enhance circulation of blood, lymphatic fluid etc. It rebalances the digestive, circulatory, respiratory, nervous and immune systems, which become unbalanced as your body redirects resources to deal with everyday stress.

In Shiatsu, the practitioner stimulates and balances the body through applying pressure and massaging acupuncture points and meridians.

This stimulates the body's self healing capability and helps to maintain or regain health.

The SMC-2000H Shiatsu massage cushion has rotating massage heads which imitates the pressure and kneading motions of a real Shiatsu massage.

3 Safety advice

The device is only intended for the purpose described in this user guide. The manufacturer cannot be held liable for personal injuries or damages occurring as a result of inappropriate or careless use.

3.1 General

- Always read the safety instructions carefully and keep this user guide for future reference.
- Follow the basic safety precautions for all electronic equipment when using this product.
- This product is not a toy. Don't let children play with it.
- The device is only for private use and is not intended for medical or commercial purposes. If you have any concern regarding your health, consult your doctor before using this product.
- This unit is not waterproof. Don't expose this product to rain or moisture.
- Do not leave the product unattended when plugged in. Unplug the product immediately from the main power after use.
- Do not place the device on delicate surfaces, as otherwise the zip fastener could potentially damage the surface.
- Make sure that the appliance is not too close to heat sources such as ovens, radiators or direct sunlight.

3.2 Electric shock

- Do not operate this product when the cord or plug is damaged.
- Do not use the product if it is not working properly, dropped or damaged.
- Repairs to electrical appliances must only be carried out by qualified persons. Incorrect repairs can lead to considerable danger for the user. In the event of repairs, please contact our customer service or an authorized dealer.
- Don't disassemble this product. It doesn't contain parts that can be repaired by the user.
- Do not connect other attachments, as the power adapter, than the one supplied with the product.
- Do not drop or insert object into any opening.
- Do not place the basic unit in a damp room and at a distance of less than 1.5 m away from a water source.
- Do not crush and avoid sharp folds.
- Do not use the product for more than 15 minutes at a time. Extensive use could lead to the product's excessive heating and shorter lifetime.
- The zip fastener on the Shiatsu massage device must not be opened. It is only fitted for reasons relating to the manufacturing of the device.
- Don't disassemble this product. It doesn't contain parts that can be repaired by the user.
- Never use the appliance during a thunderstorm.
- Do not use the massage cushion under a blanket or pillow. Excessive heating can occur and cause fire, electrocution and injuries to person.

3.3 Massage cushion precautions

- It is not allowed to use the massage cushion:
 - On young children under the age of 16.
 - On persons suffering from a herniated disk or from abnormal changes to the vertebra.
 - Not in case of cardiovascular diseases and headaches.
 - Not in the case of swelling, burns, inflammation, open wounds, fresh wounds, hematomas, bleeding or other injuries to the back or legs.
 - During pregnancy.
 - In bed or while sleeping.
 - After taking medication or consuming alcohol (reduced levels of awareness).
 - On a person with insensitive skin or on a person with poor blood circulation.
 - In a car
 - While operating machinery.
 - If you are fitted with a cardiac pacemaker, implant or other aid. Individuals with pacemaker should consult their doctor before use.
 - If you suffer from any of the following: common cold accompanied by fever, varicose veins, thrombosis, phlebitis, jaundice, diabetes, nervous system diseases (e.g. sciatica), tuberculosis, tumors, hemorrhoids, bruises or acute inflammation.
 - Not in case of pains of unclarified origin.
- Close supervision is necessary when the appliance is used by or near invalids or disabled persons.
- Do not use on this product before sleeping. The massage has a stimulating effect and can delay sleep.
- When the unit is working, do not touch the spinning massage mechanism with any other body part except your back.
- Only gentle force should be exerted against the unit in order to eliminate any/all risk of injury.
- If the material or fabric covering the massage mechanism becomes torn or damages, discontinue use.
- Never put any body part (fingers or hands) inside the interior of the unit.

4 Installation

4.1 Transportation screw

Remove the transportation screw on the backside of the product before use. Not removing the screw can cause permanent damage to the massage cushion. No guarantee is accepted in this case.

The massage cushion comes with a screw located on the backside of the unit to protect the massage mechanism during shipment.

This screw must be removed before use to allow that the massager functions properly.

Properly dispose the screw and install the included black plug into the hole that was left by the screw.

4.2 Strapping system

The massage cushion can be installed to almost any chair.

Slip the straps over the seat or chair and adjust the fastening straps to secure.

The massage won't slip or slide away.

4.3 Connecting to the main power.

Only use the power adaptor supplied. Do not use other adapters, as this may damage the device. The wall outlet for the power supply must be close and accessible.

- Plug one end of the power adaptor into the 230V outlet and the other end into the adapter jack which leads from the massage cushion.

5 Button description

I Power button

- Press the power button to turn the unit on. The corresponding LED light will illuminate.
- Press the power button once again to turn the unit off. The massage mechanism will return to its lowest point before switching off.

II Width adjustment


During rolling massage, you can change the distance between the massage heads.

With each press on the width button, the massage heads will move to the next available position.

III Full back rolling massage

Rolling is like two hands or thumbs working up and down along your spine.

- Press this button for a full back 'rolling' massage.
- The corresponding LED light will illuminate.
- You will enjoy a full back rolling massage
- Press this button again to switch this massage off.

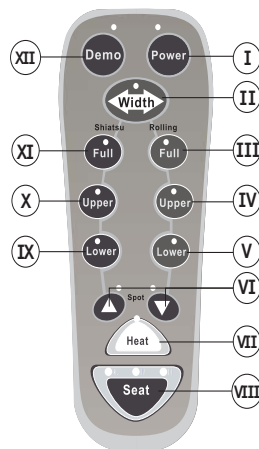IV Upper back rolling massage

- Press this button for a massage of your upper back only.
- The corresponding LED light will illuminate.
- The massage mechanism will skip the lower bag during massage.
- Press this button again to switch the massage off.

V Lower back rolling massage

- Press this button for a massage of your lower back only.
- The corresponding LED light will illuminate.
- The massage mechanism will skip the upper bag during massage.
- Press this button again to switch the massage off.

VI Spot massage (shiatsu)

- Press the spot massage up or down button to stop the kneading massage at a certain point.
- Hold the up or down button to move the massage mechanism up or down.
- Release the button when you have reached the desired massage point.

VII Heat function

- Press the heat button to activate a suitable heating function and the corresponding LED light will illuminate.
- Press the button again to turn the heat function off.

The heat function cannot be turned on if no massage program is selected.

VIII Seat vibration

- Press this button to activate the seat vibration. The vibration will turn On in low-intensity and the corresponding LED light (L) will illuminate.
- Press this button repeatedly to switch to medium (M) intensity, high (H) intensity or to deactivate the seat vibration.

IX Full back Shiatsu massage

- Press this button for a full back 'Shiatsu' massage.
- The corresponding LED light will illuminate.
- You will enjoy a full back kneading massage
- Press this button again to switch this massage off.

X Upper back Shiatsu massage

- Press this button for a massage of your upper back only.
- The corresponding LED light will illuminate.
- The massage mechanism will skip the lower bag during massage.
- Press this button again to switch the massage off.

XI Lower back Shiatsu massage

- Press this button for a massage of your lower back only.
- The corresponding LED light will illuminate.
- The massage mechanism will skip the upper bag during massage.
- Press this button again to switch the massage off.

XII Demo button

- Press this button to demonstrate the features of the massage cushion.

There is a 15-minute auto shut off on the unit for your safety.

When the massager is switched off, the massage mechanism will always return to its lowest position.

When selecting the seat vibration together with the back massage, there will be a slight decrease of the vibration strength. This adjustment is perfectly normal and should not be interpreted as a defect.

6 Technical data

Dimensions (W/H/D)	680*145*470mm
Weight	4.5 Kg
Power supply unit voltage:	
– Input:	100-240 V/AC/50-60 Hz
– Output:	DC 12V/3A

7 Cleaning

- Disconnect the device and allow it to cool down.
- Clean it with a soft cloth or slightly damp sponge.
- Never allow that water or any other liquid enters the device or accessories.
- Never use abrasive cleaners, brushes, gasoline, kerosene, glass polish or paint thinner to clean.

8 Storage

- Do not wrap the power cord around the unit before storing.
- Store the massage cushion in its original box and in a safe, dry and cool place.
- Store it out of the reach of children.
- Avoid contact with sharp edges or pointed objects as this might damage the fabric.
- You can use the transportation screw for locking the massage mechanism. Be careful to remove it from the unit when switching it on again.

9 Disposal of the device (environment)

At the end of the product lifecycle, you should not throw this product into the normal household garbage but bring the product to a collection point for the recycling of electrical and electronic equipments. The symbol on the product, user guide and/or box indicate this.

Some of the product materials can be re-used if you bring them to a recycling point. By re-using some parts or raw materials from used products you make an important contribution to the protection of the environment. Please contact your local authorities in case you need more information on the collection points in your area.

10 Topcom warranty

10.1 Warranty period

The Topcom units have a 24-month warranty period. The warranty period starts on the day the new unit is purchased.

Consumables or defects causing a negligible effect on operation or value of the equipment are not covered.

The warranty has to be proven by presentation of the original purchase receipt, on which the date of purchase and the unit-model are indicated.

10.2 Warranty handling

A faulty unit needs to be returned to a Topcom service centre including a valid purchase note.

If the unit develops a fault during the warranty period, Topcom or its officially appointed service centre will repair any defects caused by material or manufacturing faults free of charge.

Topcom will at its discretion fulfill its warranty obligations by either repairing or exchanging the faulty units or parts of the faulty units. In case of replacement, colour and model can be different from the original purchased unit.

The initial purchase date shall determine the start of the warranty period. The warranty period is not extended if the unit is exchanged or repaired by Topcom or its appointed service centres.

10.3 Warranty exclusions

Damage or defects caused by incorrect treatment or operation and damage resulting from use of non-original parts or accessories not recommended by Topcom are not covered by the warranty.

The warranty does not cover damage caused by outside factors, such as lightning, water and fire, nor any damage caused during transportation.

No warranty can be claimed if the serial number on the units has been changed, removed or rendered illegible.

**This product is in compliance with the essential requirements and other relevant provisions of the directive 2004/108/EC.
The Declaration of conformity can be found on:
<http://www.topcom.net/cedeclarations.php>**

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