# BébéSounds®

# Digital Baby Scale

Balanza Digital para Bebé

Model #BR149

**Instruction Manual**  
Manual de Instrucciones

**UNISAR®**

[www.bebesounds.com](http://www.bebesounds.com)

## TABLE OF CONTENTS

|                                      |     |
|--------------------------------------|-----|
| Introduction.....                    | 4   |
| Battery Installation.....            | 5   |
| Instructions for Use.....            | 5   |
| Care and Cleaning.....               | 6   |
| Technical Specifications.....        | 6   |
| Limited Product Warranty.....        | 7   |
| Growth Charts.....                   | 8-9 |
| Also Available from Bébésounds®..... | 10  |
| En Español.....                      | 11  |


## INTRODUCTION

Thank you for purchasing our Bébésounds® Digital Baby Scale. As you know, it is important to monitor the progress of your baby's physical development. Our Digital Baby Scale will enable you to track your baby's growth by accurately measuring her weight and length.

We have included two growth charts published by the National Center for Health Statistics—one for girls and one for boys. After weighing your baby you can use these charts to determine what percentile she or he is in based on the national average. For example, if your baby is in the 90th percentile for weight, that means 90% of babies of the same sex and age group in the United States weigh the same or less than she does and 10% weigh more. However, please keep in mind that these are only averages and should only be used as a general guide. Each baby develops at her own individual pace. Please consult your physician if you have any concerns.

### Key Features


- Weight is measured in pounds, ounces, or kilograms
- Measures baby's length in inches and centimeters
- Jumbo, easy-to-read LCD display
- Memory stores and recalls last weight measurement
- Special Tare feature deducts weight of blanket if used when weighing your baby
- Records weight of up to 44 pounds/704 ounces/20 kilograms
- Accurately measures weight in increments of 0.02 pounds/ 0.5 ounces/0.01 kilograms
- Measures length of up to 20.5 inches/52 centimeters
- Low battery indicator
- Automatic shut-off ensures long battery life
- Smooth, contoured platform is comfortable for baby
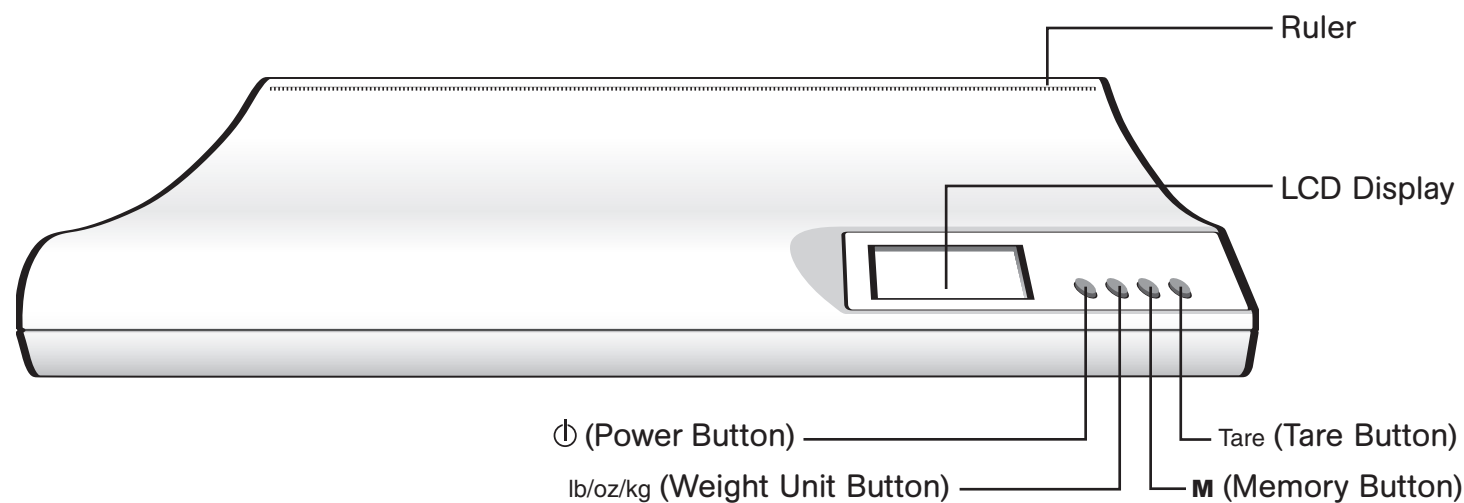
## BATTERY INSTALLATION (requires 4 AAA alkaline batteries, included)

Use only AAA alkaline batteries. Always replace all batteries at the same time. Do not mix new and old batteries. For long periods of non-use, remove batteries from the compartment.




1. Remove the battery compartment cover on the bottom of the scale.
2. Insert four (4) AAA alkaline batteries matching positive (+) to positive (+) and negative (-) to negative (-).
3. Replace the battery cover.
4. Replace the batteries when the low battery indicator appears (see Fig. 1).
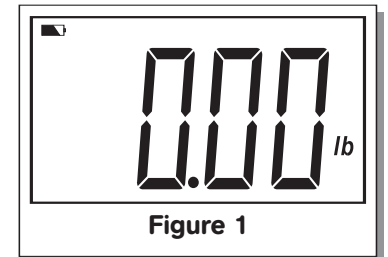
## INSTRUCTIONS FOR USE

### Without a Blanket


1. Lay the scale on a flat, firm surface.
2. Turn the scale on by pressing the Power button (⏻).
3. The last weight measurement saved in the memory is displayed in pounds. To see this weight in ounces or kilograms, press the Weight Unit (lb/oz/kg) button.
4. After a few seconds, the LCD will display 0.00 and "Ready" will appear (see Fig. 2). Press the Weight Unit (lb/oz/kg) button if a different measurement unit is desired.
5. Put your baby on the scale. When "Hold" flashes, the weight measurement is complete (see Fig. 3).
6. To store the measurement into memory, press the Memory button (M). "M" will then flash and the measurement is stored (see Fig. 4).
7. To take an additional measurement, turn the scale off and then on again. Repeat steps above.
8. Press the Power button (⏻) to turn the scale off or wait approximately 10 seconds and it will turn off automatically.

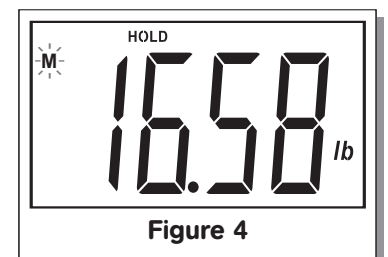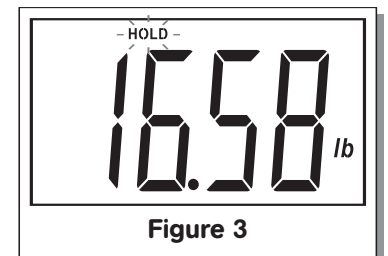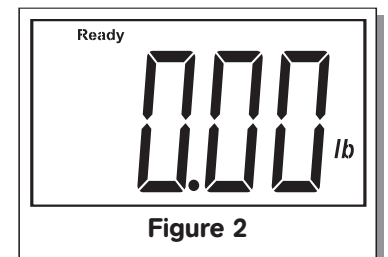
### With a Blanket

Our special Tare feature allows you to weigh your baby with a blanket if desired. It will deduct the extra weight (if any) from the blanket. To weigh your baby with a blanket or any other additional item, please follow these instructions:

1. Follow steps 1 through 4 as listed above.
2. Put the blanket on the scale. When "Hold" flashes, the weight measurement is complete. Press the Tare button.
3. The scale will display 0.00 and "Ready" and "Tare" will appear (see Fig. 5). If you remove the blanket from the scale, a negative weight will be displayed (see Fig. 6).

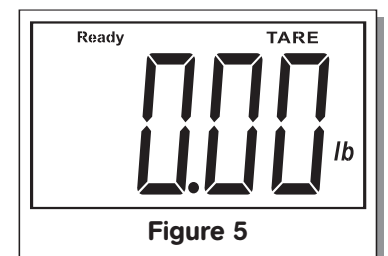4. Put your baby on the scale with the blanket. When “Hold” flashes, the weight measurement is complete.
5. To store the measurement into memory, press the Memory button (M). “M” will then flash and the measurement is stored.
6. To take an additional measurement, turn the scale off and then on again. Repeat steps above.
7. Press the Power button (⏻) to turn the scale off or wait approximately 10 seconds and it will turn off automatically.

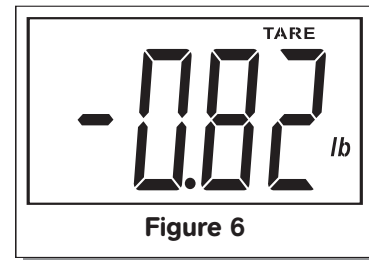**WARNING:** Never leave your baby unattended during the weighing process.

## CARE AND CLEANING

1. To clean the scale, simply wipe with a soft dry or slightly damp cloth. Only use a mild soap, never abrasives.
2. Do not submerge in water or any other liquid.
3. Store in a cool dry place.
4. Keep out of reach of children when not in use.
5. During long periods of non-use, remove batteries from the compartment.

## TECHNICAL SPECIFICATIONS

- For home use only
- Requires four (4) AAA alkaline batteries (included)
- Maximum weight capacity: 44 pounds/704 ounces/20 kilograms
- Maximum length capacity: 20.5 inches/52 centimeters
- Graduation: 0.02 pounds/0.5 ounces/0.01 kilograms

## APPLIED STANDARDS AND APPROVALS

This product complies with the requirements of the EC Council Directive on electromagnetic compatibility 89/336/EEC.

The following standards apply to design and/or manufacturing of the product:

- EN61000-6-3:2001: Electromagnetic compatibility. Generic emission standard - Part 1: Residential, commercial and light industry
- EN61000-6-1:2001: Electromagnetic compatibility. Generic immunity standard - Part 1: Residential, commercial and light industry
- EN61000-3-2:2000: Electromagnetic compatibility, Part 3 limit. Section 2, limits for harmonic current emissions (equipment input current  $\leq 16$  A per phase)
- EN61000-3-3:1995: Electromagnetic compatibility. Limits of voltage changes and fluctuations and flicker in low-voltage supply systems for equipment with rated current  $\leq 16$  A per phase

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two condi-

tions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

## LIMITED PRODUCT WARRANTY

Register your product warranty online at [www.bebesounds.com](http://www.bebesounds.com).

Unisar Inc., the distributor of the Bébésounds® Digital Baby Scale, warrants to the original owner of this product that this product is free from defects in materials and workmanship for a period of one (1) year from the date of the original purchase. If this product fails to function properly when used as directed under normal conditions within the one-year period, Unisar Inc. will either repair or replace the product, at its option, free of charge, when it is returned in its original packaging, postage prepaid, with proof of purchase. Unisar Inc. will bear the cost of repairing or replacing the product and shipping it back to you within the United States. Repair or replacement of the defective item is subject to verification of the malfunction or defect when delivered to:

Unisar® Service Center  
15 West 36th Street  
New York, NY 10018


NOTE: The warrantor expressly disclaims liability for incidental, consequential, or special damages of any nature. This warranty does not cover 1) A defect that has been caused by repair or modification by any person other than a Unisar® technician; 2) A defect caused by negligence, accident, or misuse of the product; 3) A defect caused by incorrect installation of and/or defective batteries; 4) Battery life; or 5) Damage resulting from improper operation.

LEGAL RIGHTS VARY FROM STATE TO STATE.

NOTE: This warranty is valid only for products sold in the United States. If you have purchased your product outside of the United States, please see your local retailer or distributor for service.

If you have any questions, please call our Help Desk toll free at 1-888-232-6476 Mon-Fri from 9:00AM-5:00PM EST.

## CDC Growth Charts: United States

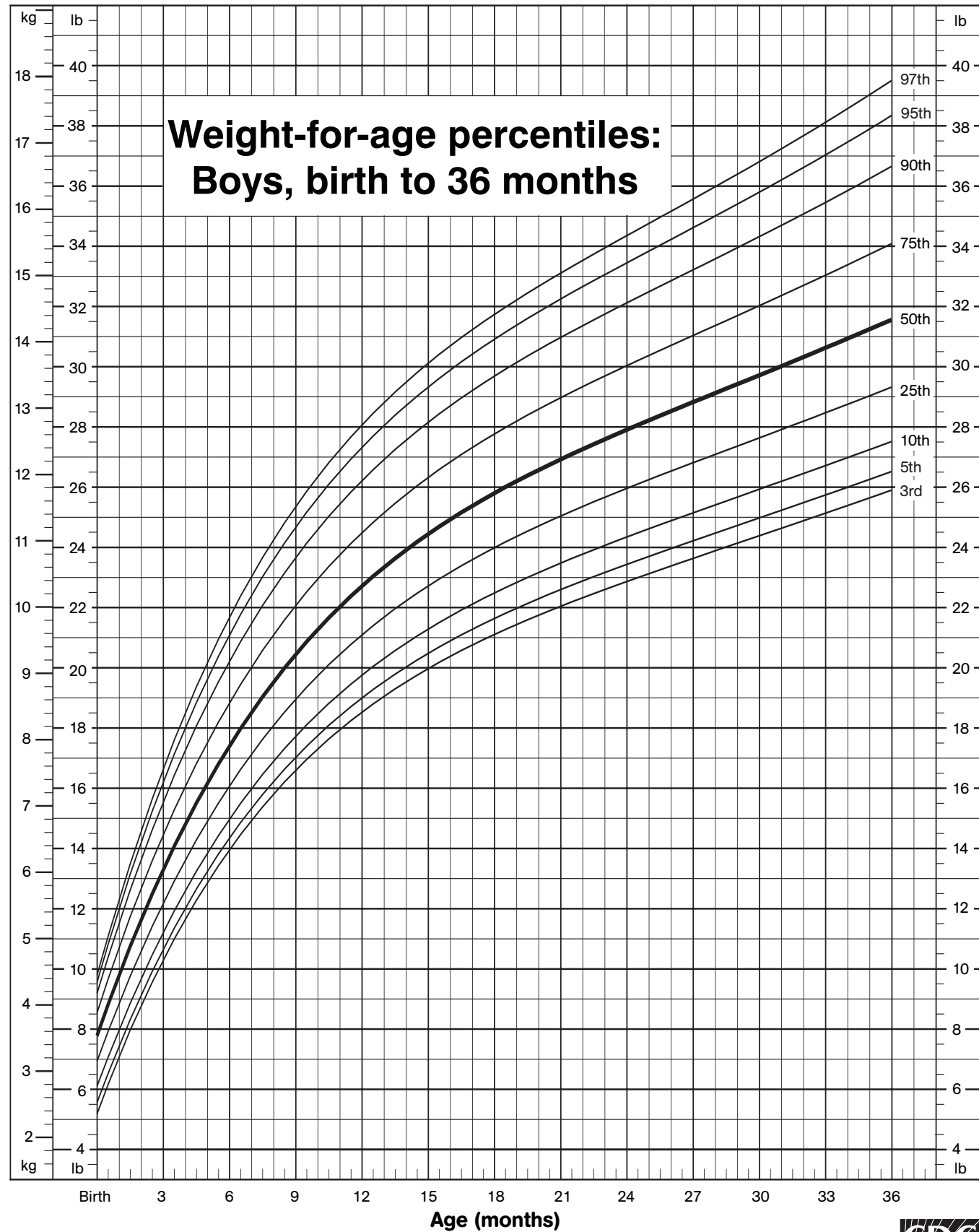
Published May 30, 2000.  
SOURCE: Developed by the National Center for Health Statistics in collaboration with the National Center for Chronic Disease Prevention and Health Promotion (2000).

SAFER • HEALTHIER • PEOPLE™

## CDC Growth Charts: United States

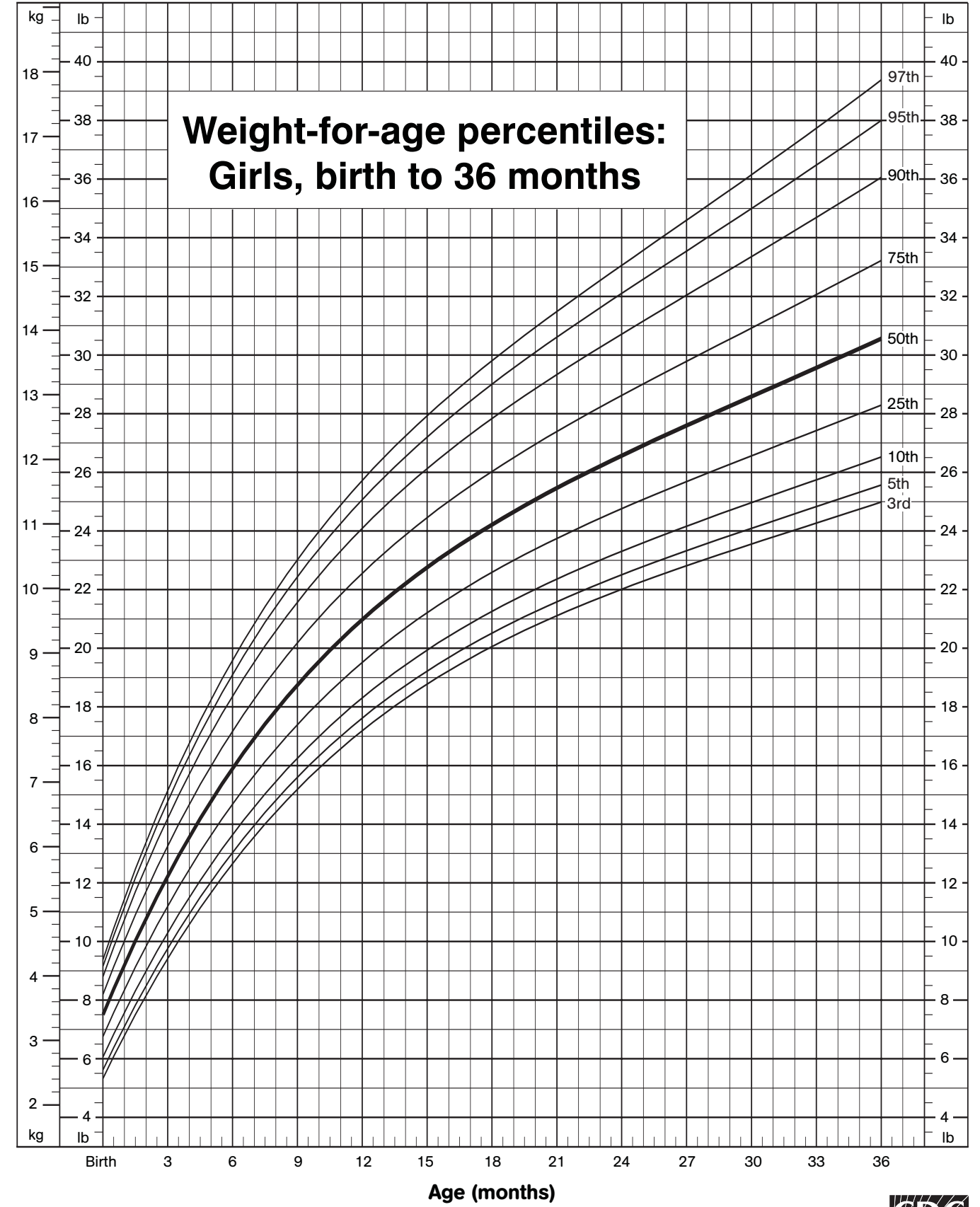Published May 30, 2000.  
SOURCE: Developed by the National Center for Health Statistics in collaboration with the National Center for Chronic Disease Prevention and Health Promotion (2000).

SAFER • HEALTHIER • PEOPLE™

**NasalClear™ Aspirator**

Safely and quickly clear your baby's nasal congestion with this battery operated aspirator. Gentle enough to use on a newborn baby. Comes apart for easy and thorough cleaning.

**AlwaysClean™ Pacifier**

Closes when dropped! Unique shield automatically closes when the pacifier is dropped keeping it sanitary for your baby. Silicone orthodontic nipple. 0-6 months and 6+ months.

**Angelcare® Movement Sensor with Sound Monitor**

Stop worrying while your baby sleeps!™ Detects your baby's slightest movements, even while she sleeps. An alarm sounds ONLY if absolutely no movement is detected for 20 seconds. Also works as a regular sound monitor.

**1-Second Ear Thermometer**

Makes taking your child's temperature quick and convenient. No probe cover or "ear-tug" required. Cleans with alcohol. Lighted jumbo display. Measures in Fahrenheit and Celsius

**5-Second 3-Way Flexible Thermometer**

Get an accurate temperature reading in just 5 seconds! And the flexible tip provides your child with maximum comfort. Recalls last temperature reading. Measures in Fahrenheit and Celsius.

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©2006 Listro Associates