
# Oral-B®

## AdvancePower

### D 9513 K

## Kids

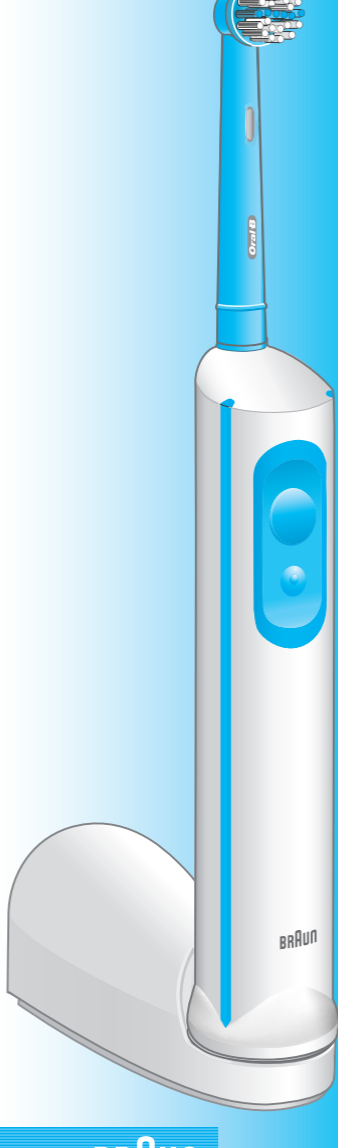powered by **BRAUN**

Type 4728

English

Thank you for purchasing an Oral-B AdvancePower Kids, which is designed and engineered by Braun.

If you have any questions, please call:

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





Printed in Germany

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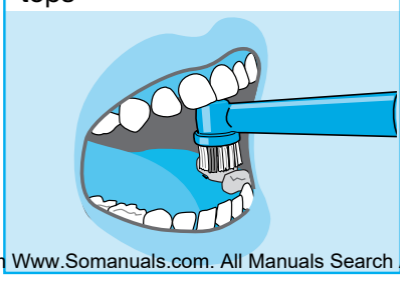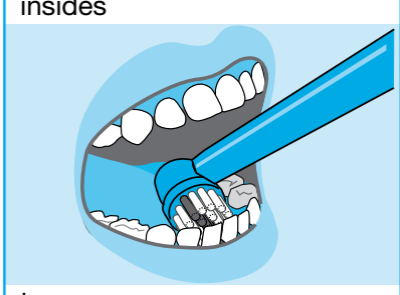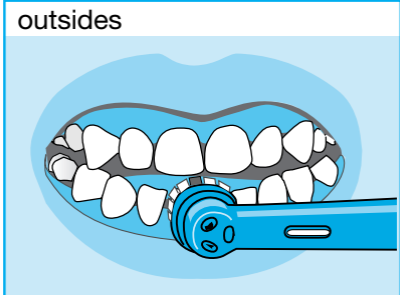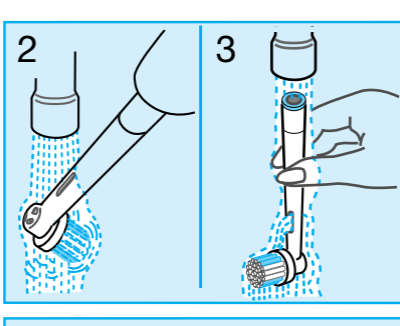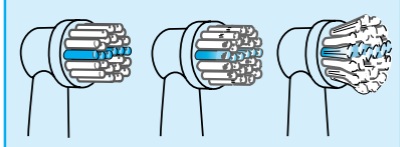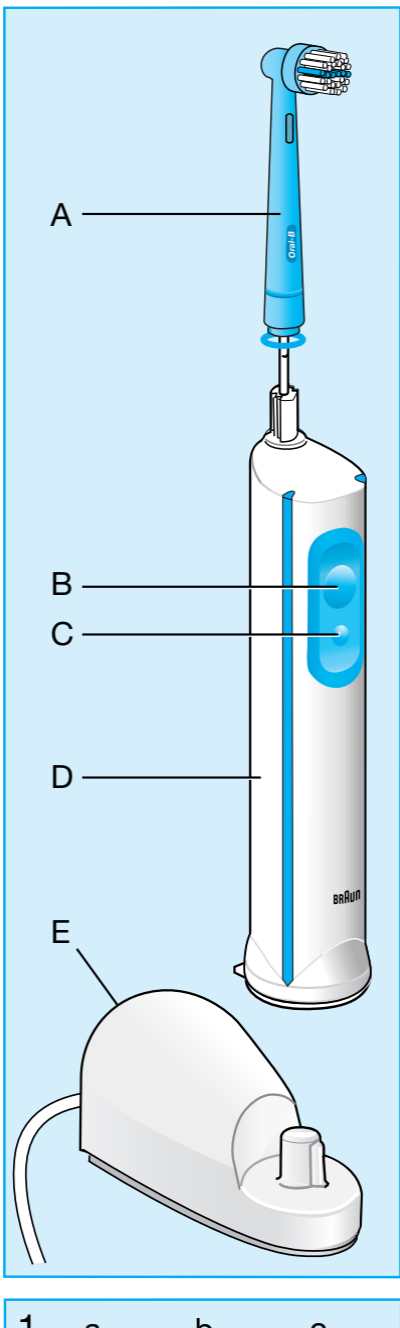# IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

## READ ALL INSTRUCTIONS BEFORE USING.

## DANGER

To reduce the risk of electrocution:

1. Do not place or store charger where it can fall or be pulled into a tub or sink.
2. Do not place charger in or drop into water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

---

## WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons:

1. Close supervision is necessary when this product is used by, on, or near children.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charger has been dropped into water. Return the product to a service center for examination and repair.
4. Keep the cord away from heated surfaces.
5. Never drop or insert any object into any opening.

## SAVE THESE INSTRUCTIONS

The Oral-B AdvancePower Kids has been developed to effectively remove plaque from children's teeth. We hope your child will enjoy his/her new Oral-B appliance.

This appliance is suitable for children 3 years or older. Please read the use instructions carefully and completely before showing your children how to use the appliance. We recommend that parents supervise their children when they start using the appliance in order to help them to become familiar with it.

#### Important

Periodically check the cord for damage. Should this be the case, take the charging unit (E) to a Braun Service Center. A damaged or non-functioning unit should no longer be used.

### Description

- A Brushhead
- B On switch
- C Off switch
- D Handle
- E Charging unit

### Specifications

Power supply: see bottom of charging unit

Handle voltage: 1.2 V

### Connecting and charging

The toothbrush is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit. The handle must be switched off.
- A full charge takes 16 hours and will provide an operation time of approximately 20 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

#### Battery maintenance

To maintain the maximum capacity of the rechargeable battery, you must discharge it at least once every 6 months. To do this, unplug the charging unit and allow the battery to run down with normal use.

### Using the toothbrush

Use the following steps to show and explain to your child how to use the appliance properly:

- Put a small amount of toothpaste onto the brushhead.
- To avoid splashing, guide the brushhead to the teeth, before switching the appliance on.
- Start with the lower teeth: Slowly guide the brushhead from tooth to tooth;
  - first the outsides,
  - then the insides and
  - finally the chewing surfaces.
- The Oral-B AdvancePower Kids automatically provides the proper cleaning action. Do not press too hard, nor scrub, simply let the brush do all the work.
- After one minute of operation, an incorporated timer signals by a short stutter. This tells your child to move on the upper teeth and clean them – again tooth by tooth – until the next timer signal is given. This way you can make sure that your child brushes his/her teeth for at least 2 minutes as recommended by dental professionals.
- Turn off the handle before taking the brushhead out of the mouth.

Dentists recommend that you replace any toothbrush every 3 months. The brushhead is provided with a horizontal row of blue INDICATOR® bristles (Ia). With proper brushing twice a day for two minutes, the color will fade half-way down the bristles in a 3-month period (Ib).

If bristles begin to splay before the color recedes, this indicates that your child exerts too much pressure on teeth and gums (Ic).

#### Replacement brushheads

EB10 replacement brushheads for children are available at your retailer or at Braun Service Centers.

### Cleaning

After use, rinse the brushhead for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brushhead. Rinse both parts separately under running water (III) and wipe them dry.

### Environmental notice

This product contains nickel-cadmium rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Center.

Ni-Cd

However, if you prefer to dispose of the batteries yourself, please refer to the section «Battery removal».

Subject to change without notice

from [www.Somanuals.com](http://www.Somanuals.com). All Manuals Search And

---

# FOR USA ONLY

---

## Limited 2-Year Warranty

Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

## To obtain service in the USA

Should this product require service, please call 1-800-BRAUN 11 (272-8611) to be referred to the authorized Braun Service Center closest to you.

*To obtain service:*

- A Take the product in to the authorized Braun Service Center of your choice or,
- B Ship the product to the authorized Braun Service Center of your choice.
  - Pack the product well.
  - Ship the product prepaid. Insurance is recommended.
  - Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts and accessories:

Contact an authorized Braun Service Center. Please call 1-800-BRAUN 11 (272-8611) to be referred to the authorized Braun Service Center closest to you.

To order original Braun replacement parts and accessories (including replacement brushheads) by phone:

Please call 1-800-892-8534.

To speak to a Braun Consumer Service Representative:

Please call 1-800-984-9668.

---

# FOR CANADA ONLY

---

## Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service centre, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

## How to obtain service in Canada

Should this product require service, please call 1 800 387-6657 to be referred to the authorized Braun Service Center closest to you.

*To obtain service:*

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice.  
Pack the product well.  
Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

*To obtain original Braun replacement parts & accessories:*

Contact an authorized Braun Service Center. Please call 1 800 387-6657 to be referred to the authorized Braun Service Center closest to you.

*To speak to a Braun Service Representative:*

Please call 1 800 387-6657.

---

# 12 Month Performance Guarantee

---

## For USA only

Use the Oral-B AdvancePower Kids daily until you are fully familiar with this mechanical plaque removal instrument. If you are not convinced that the Oral-B AdvancePower Kids significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Oral-B AdvancePower Kids within 12 months, post-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

### 12 Month Performance Guarantee

Braun Division,  
The Gillette Company,  
Sobin Park and A Street, MS 1F  
South Boston, MA 02127

**Sorry:** No refunds will be made for any appliance damaged by accident, neglect or unreasonable use.

The Oral-B AdvancePower Kids must be shipped in its original packaging and be postmarked no later than 12 months after the date of purchase.

Please allow 6–8 weeks for handling.

---

Name of Regular Toothbrush

---

Your Name

---

Street

---

City

---

State

Zip Code

---

# 12 Month Performance Guarantee

---

## For Canada only

Use the Oral-B AdvancePower Kids daily until you are fully familiar with this mechanical plaque removal instrument. If you are not convinced that the Oral-B AdvancePower Kids significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Oral-B AdvancePower Kids within 12 months, postage-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

### 12 Month Performance Guarantee

Braun Products  
Gillette Canada Company  
4 Robert Speck Parkway  
Mississauga ON L4Z 4C5

**Sorry:** No refunds will be made for any appliance damaged by accident, neglect or unreasonable use.

The Oral-B AdvancePower Kids must be shipped in its original packaging and be postmarked no later than 12 months after the date of purchase.

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## Nickel-Cadmium Battery Pack

### SAVE THIS INFORMATION

PRODUCT CONTAINS A NICKEL-CADMIUM RECHARGEABLE BATTERY PACK. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

May not be disposed of in municipal waste. Recycle or dispose of in an environmentally sound manner.

**RBRC**

Ni-Cd

**Ni-Cd**

Call Braun 800-272-8611 for recycling information.

See next page for easy battery pack removal.

DISPOSE OF YOUR BATTERY PACK PROPERLY.

DO NOT PLACE USED BATTERIES IN YOUR REGULAR TRASH.

THIS NICKEL-CADMIUM BATTERY PACK MUST BE COLLECTED, RECYCLED OR DISPOSED OF IN AN ENVIRONMENTALLY SOUND MANNER.

The incineration, landfilling or mixing of nickel-cadmium batteries with the municipal solid waste stream is PROHIBITED BY LAW in most areas.

Return this battery pack to a federal or state approved nickel-cadmium battery recycler. This may be where you purchased the battery.

In MINNESOTA call 1-800-225-PRBA if further disposal information is required; or Contact your local waste management officials for other information regarding the environmentally sound collection, recycling and disposal of this battery pack.

For more information, visit [www.braun.com](http://www.braun.com).

© 2011 Braun. All rights reserved. Printed in the USA.

See next page for easy battery pack removal.

DISPOSE OF YOUR BATTERY PACK PROPERLY.

DO NOT PLACE USED BATTERIES IN YOUR REGULAR TRASH.

THIS NICKEL-CADMIUM BATTERY PACK MUST BE COLLECTED, RECYCLED OR DISPOSED OF IN AN ENVIRONMENTALLY SOUND MANNER.

The incineration, landfilling or mixing of nickel-cadmium batteries with the municipal solid waste stream is PROHIBITED BY LAW in most areas.

Return this battery pack to a federal or state approved nickel-cadmium battery recycler. This may be where you purchased the battery.

In MINNESOTA call 1-800-225-PRBA if further disposal information is required; or Contact your local waste management officials for other information regarding the environmentally sound collection, recycling and disposal of this battery pack.

For more information, visit [www.braun.com](http://www.braun.com).

© 2011 Braun. All rights reserved. Printed in the USA.

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Gillette Distribuidora, S.A. de C.V.