Oral-B®

Professional Care Center

**OC 17545X
OC 17525**

powered by **BRAUN**

Type 3719

Deutsch	4, 103	(D) (A)	00 800 27 28 64 63
English	9, 103	(CH)	08 44 - 88 40 10
Français	14, 103	(GB)	0 800 73 11 792
Español	20, 103	(IRL)	1 800 509 448
Português	25, 103	(F)	0 810 365 855
Italiano	30, 104	(B)	0 800 14 592
Nederlands	35, 104	(LU)	4 00 50 51
Dansk	40, 104	(E)	901 11 61 84
Norsk	45, 104	(P)	808 20 00 33
Svenska	50, 105	(I)	02-667 86 23
Suomi	55, 105	(NL)	0800-4 45 53 88
Polski	60, 105	(DK)	70 15 00 13
Ceský	65, 106	(N)	22-72-88-10
Slovenský	70, 106	(S)	020 - 21 33 21
Magyar	75, 106	(FIN)	020 377 877
Ελληνικά	80, 106	(PL)	0 801 127 286 0 801 1 BRAUN
Türkçe	86,	(CZ) (SK)	00420-800-11-33-22
Русский	91, 107	(H)	+36 (1) 801-3800
عربي	102, 96	(GR)	210 9478700
		(TR)	0 212 473 75 85
		(RUS)	+7 495 258 62 70

Internet:

www.oralb.com
www.braun.com

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The Oral-B Professional Care Center has been developed with leading dental professionals to effectively remove plaque from all areas of your teeth and to attack plaque bacteria.

The Oral-B Professional Care toothbrush adds a gentle pulsating movement to its oscillation to create a deep-cleaning 3D brushing action. The toothbrush should be used twice a day for at least two minutes, the minimum brushing time recommended by dental professionals.

The Oral-B OxyJet® irrigator mixes air into water in the form of micro-bubbles. These micro-bubbles in the water stream attack plaque bacteria. At the same time, the irrigator removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

The appliance must be mounted in a way that it cannot fall into water.

The plug should always be accessible. Periodically check the cord for damage. A damaged or non-functioning appliance should no longer be used. Take it to a Braun Service Centre.

Description

Toothbrush

- ① Charging base
- ② Charge indicator
- ③ Toothbrush handle (rechargeable)
- ④ On/off switch
- ⑤ Brushhead
- ⑤a Interspace brushhead (only with OC 17545X model)
- ⑥ Brushhead compartment

Irrigator

- ⑦ Jet compartment
- ⑧ Irrigator jet
- ⑨ Jet switch
- ⑩ Jet release button
- ⑪ On/off switch
- ⑫ Water pressure control dial
- ⑬ Handle switch (water stop and go)
- ⑭ Irrigator handle
- ⑮ Water container (600 ml)

Specifications

Power supply: see bottom of housing

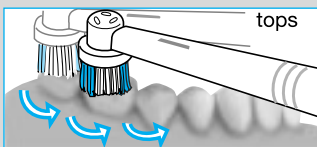
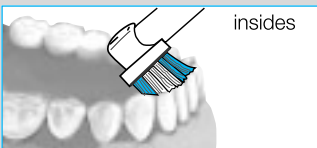
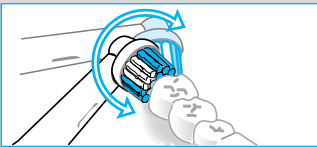
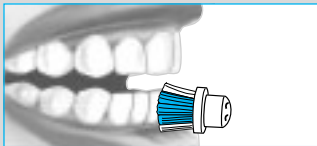
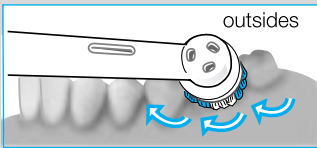
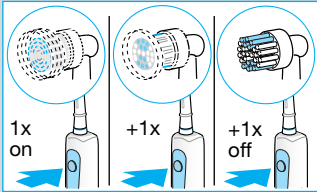
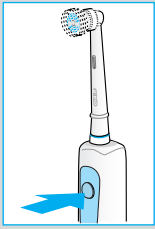
Power intake: 24 W

Using the toothbrush

Connecting and charging

The toothbrush is electrically safe and can be used in the bathroom without hesitation.

- Plug the appliance into an electrical outlet. Place the toothbrush handle ③ on the charging base ①.
- The charge indicator ② flashes while the handle is being charged.
- When the handle is fully charged, the charge indicator lights permanently. A full charge takes at least 16 hours and provides an operation time of approximately 45 minutes.
- For everyday use, the toothbrush handle can be stored on the charging base to maintain it at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, the following procedure should be done every 6 months: Do not place the toothbrush handle back onto the charging base in order to discharge it by regular use. Once discharged, recharge the handle on the charging base.

Operating the on/off switch ④

The toothbrush features two speed options, the lower speed being ideal for sensitive areas in the mouth.

- Press the switch ④ to operate with normal speed,
- press a second time to operate with reduced speed,
- press a third time to switch the toothbrush off.

Removing plaque

The toothbrush can be used with any leading brand of toothpaste.

- To avoid splashing, guide the brushhead to your teeth before switching on.
- Guide the brushhead slowly from tooth to tooth. Hold the brushhead in place for a few seconds before moving on the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces. Do not press too hard nor scrub, simply let the brush do all the work.

During the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or hygienist.

Pressure sensor

For safety and optimum brushing results, the toothbrush has a pressure control system: If too much pressure is applied, the oscillating movement of the brushhead will continue, but the pulsation will stop. There is a noticeable

change in feeling and sound when the pressure sensor is activated.

Professional timer

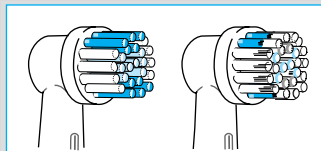
With a signal you can hear and feel, the timer invites you to brush quadrant by quadrant. After every 30 seconds, it signals by a very short stutter that it is time to move to the next quadrant. After two minutes – which is the recommended brushing time – the timer signals by a longer stutter.

INDICATOR® bristles

The INDICATOR® bristles help you monitor and improve the brushing performance.


With proper brushing, twice a day for two minutes, the blue colour on the INDICATOR® bristles will disappear halfway down within 3-4 months, signalling it is time to replace the brushhead.

Interspace brushhead

(available with OC 17 545 X model)

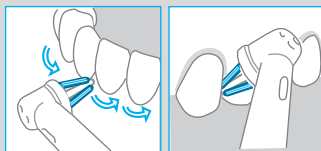The interspace brushhead is designed for cleaning between teeth as well as around bridges, crowns and implants. To use, position brushhead at 45° angle to gumline. Slowly guide the brushhead along the gumline, pausing between teeth for optimum results.


Replacement brushheads

Available at your retail store or Braun Service Centres.

After use

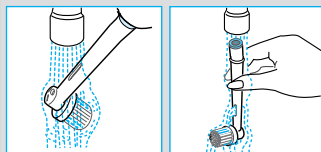
Rinse the brushhead thoroughly under running water for several seconds with the handle switched on. Switch off the handle and remove the brushhead. Clean both parts separately under running water, then wipe them dry.

Environmental notice



The toothbrush handle is provided with a nickel-hydride rechargeable battery, which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Centre.

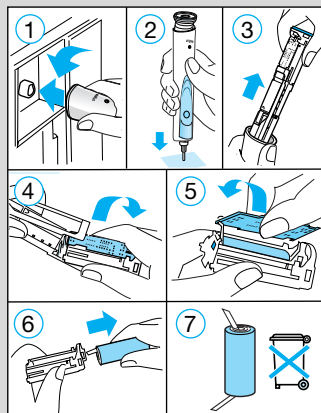However, if you prefer to dispose of the battery yourself, open the handle as shown, remove the battery and dispose of it according to local environmental regulations. Caution: Opening the handle will definitely destroy it.



Using the OxyJet® irrigator

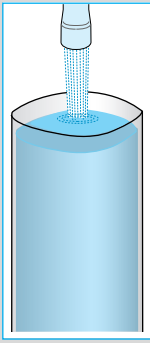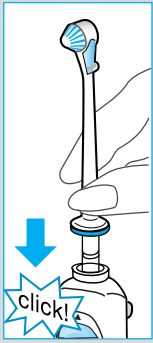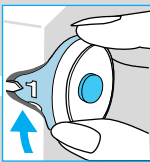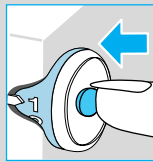Before using for the first time, fill the water container and allow the water to run through completely by holding the irrigator handle over a washbasin.

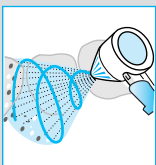
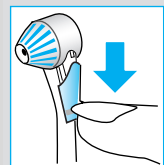
1. Attach one of the irrigator jets to the handle until it locks into place. Fill the water container with water and add mouthwash, if desired.

2. Turn on the irrigator by pressing the on/off switch. (When doing so, the handle switch should be in the lower position = water stop).
3. Select the water pressure with the pressure control dial (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.

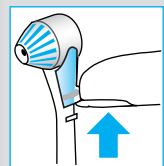
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch  up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.

You can choose between two different jet functions by sliding the jet switch  down or up (before doing so, turn off the water stream by sliding the handle switch  down):

Rotating stream for general cleaning and gum massage (jet switch in lower position): The rotation of the water stream and its micro-bubbles is designed for easy penetration below the gumline.

Mono stream: straight water stream for site-specific cleaning (jet switch in upper position)

Max. operating time: 15 minutes
Cooling off time: 2 hours

To ensure the jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

After use

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended the use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator. It will automatically turn off after 10 minutes if you forget to turn it off manually. In this case, to turn it on again, you have to press the on/off switch twice.

To take off the irrigator jet, press the jet release button. The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

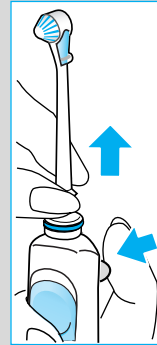
Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives
EMC 89/336/EEC and Low Voltage 2006/95/EC.

At the end of the product's useful life, please dispose of it at appropriate collection points provided in your country.

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English

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brushheads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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Üretici firma ve CE işareti uygunluk de©erlendirme kuruluşu: