
PULSONIC

S 15.523.2

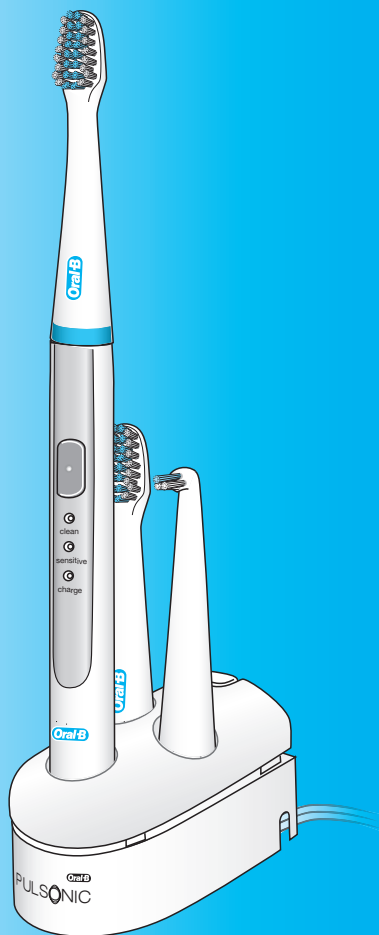Type 3715
Type 3714

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English

IMPORTANT SAFEGUARDS

Oral-B Pulsonic has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. However, when using electrical products, basic safety precautions should always be followed.

Read all instructions before using.

IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral B Braun Service Center. A damaged or non-functioning unit should no longer be used.
- This product is not intended for use by children under age three.
- Close supervision is necessary when this product is either used by or near children or invalids.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not modify the product. Do not disassemble or repair. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service center. Do not disassemble the product except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Never drop the appliance or insert any object into any opening of the appliance.
- Do not touch the power plug with wet hands. This can cause electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments or chargers not recommended by the manufacturer.
- Keep cord away from heated surfaces.
- Do not put metal objects such as coins and clips on the charger. This may generate heat damaging the unit and heating the metal item presenting a potential burn hazard.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

- a Brush head
- b Precision Tip
- c On/off button
- d Mode indicators
- e Charge indicator
- f Handle
- g Charging unit with handle stand and release button

Connecting and charging

The toothbrush is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit (g) into an electrical outlet and place the handle (f) on the charging unit ①.
- A full charge takes at least 12 hours and will provide an operation time of approximately 30 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle by regular use at least every 6 months.

Using your toothbrush

Attaching the brush head

Slide the brush head (a) on the handle and turn until it clicks into place as shown in picture ②.

Brushing modes

Oral-B Pulsonic offers different brushing modes for varying oral care needs:

«Clean» – Gentle, effective mouth cleaning.

«Sensitive» – Gentle cleaning for sensitive areas.

Your toothbrush automatically starts in the «Clean» mode when pressing the on/off button (c) ③. To shift from the «Clean» to the «Sensitive» mode, press the on/off button once again. To turn your toothbrush off, press the on/off button.

Brushing technique

- Wet brush head and apply toothpaste. You can use it with any toothpaste.
- To avoid splashing place the brush head on your teeth, then switch on.
- Move the brush head in a slightly circular motion, similar to a manual toothbrush ④.
- After a few seconds, guide the bristles to the next section
- During the first few days of use, it is possible that your gums may bleed slightly. Should bleeding persist after 2 weeks, contact your dentist.
- Replace brush head every 3 months or when bristles are worn.

Timer

A short stuttering sound at 30-second intervals (one break), which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a better overall clean.

A long stuttering sound (two breaks) indicates when the professionally recommended 2-minute brushing time has been reached.

Precision Tip

The Precision Tip (b) is designed for cleaning between teeth as well as around bridges, crowns and implants. To use, position at a 45 degree angle to gumline and slowly guide it along the gumline from one interdental space to the next ⑤.

Cleaning

Switch off the handle and remove the brush head. Rinse both the brush head and handle under running water and wipe them dry ⑥.

From time to time also clean the charging unit with a damp cloth. For thorough cleaning you may remove the handle stand by pushing the release button on top of the charging unit.

Subject to change without notice.

This product conforms to the European Directives EMC 2004/108/EC and Low Voltage 2006/95/EC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

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Warranty (Hong Kong)

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover, damage due to improper use, normal wear or use, especially regarding the brushheads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

To obtain service within guarantee period, hand in or send the complete appliance with your sales receipt to the following authorized Oral-B Braun Customer Service Centre.

Audio Supplies Company Limited

Room 506-8 St. George's Building, 2 Ice House Street, Central,

Hong Kong (852) 2524-9377

For any enquiries, please call our Consumer Relations Hotline:

Procter & Gamble Hong Kong Limited 6/F, Shui On Centre, Wan Chai,

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