Oral-B®

Vitality™
Sonic

S 12513

S 12013

powered by **BRAUN**

Type 3709

GB	0 800 73 11 792
IRL	1 800 509 448
F	0 810 365 855
B	0 800 14 592
LU	4 00 50 51
PL	0 801 127 286 0 801 1 BRAUN
CZ	221 804 335
SK	02/5710 1135
H	06-1-451-1256
HR	01/66 90 330
SLO	080 2822
TR	0 212 473 75 85
RUS	8-800-200-20-20
UA	+38 044 428 65 05
HK	852-25249377 (Audio Supplies Company Ltd.)

Internet:

www.oralb.com

www.braun.com

www.service.braun.com

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




om Www.Somanuals.com. All Manuals Search An

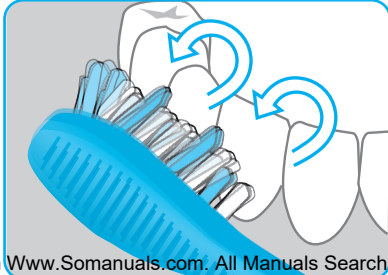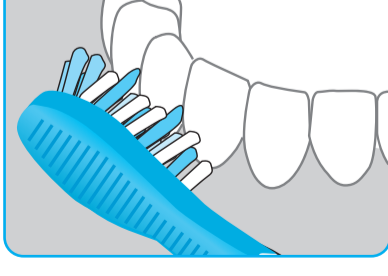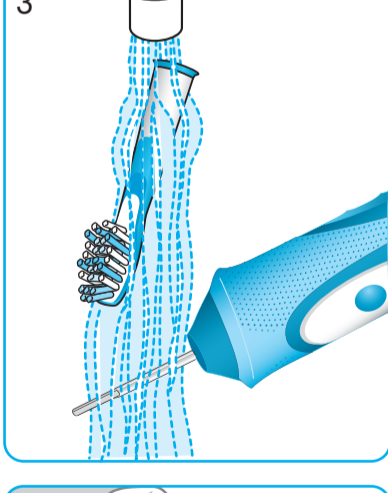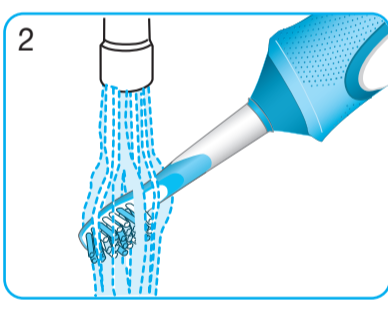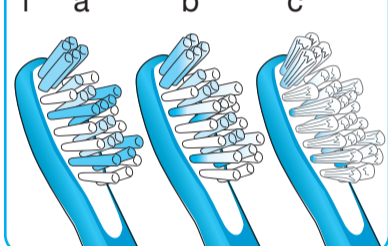English

Please read the use instructions carefully

before first use.

Important

- Periodically check the cord for damage. A damaged or non-functioning unit shall no longer be used. Take the charging unit (D) to an Oral-B Braun Service Centre.

- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are given supervision by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.

Description

A Brushhead

B On/off button

C Handle

D Charging unit

Charging

- Plug the charging unit (D) into an electrical outlet and place the handle on the charging unit.
- A full charge takes 16 hours and will provide an operation time of approximately 20 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, you must discharge it at least once every 6 months. To do this, unplug the charging unit and allow the battery to run down with normal use.

Using the toothbrush

- Slide brushhead on handle.
- Wet brushhead and apply toothpaste.
- Place brushhead on teeth, then switch on.
- Move the brushhead in a slightly circular motion, similar to a manual toothbrush. Allow the longer Interdental bristles to reach between your teeth, while the CrissCross bristles focus on cleaning the tooth surface and gumline regions.
- After a few seconds, guide the bristles to the next section.
- After two minutes (brushing time recommended by dentists) the timer (S 12513 model only) will signal by a short stutter.
- During the first few days of use, it is possible that your gums may bleed slightly. Should bleeding persist after 2 weeks, contact your dentist.
- Replace brushhead every 3 months, when bristles are worn (1c) or when colour on Indicator bristles (1b) has faded half-way down the bristles.

Cleaning

After use, rinse the brushhead under running water (2). Remove the brushhead and rinse brushhead and handle separately (3), then dry.

Subject to change without notice.

This product conforms to the European

Directives EMC 2004/108/EC and Low

Voltage 2006/95/EC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose

of the product in the household waste

at the end of its useful life. Disposal can take

place at an Oral-B Braun Service Centre or at

appropriate collection points provided in your

country.

Guarantee

We grant 2 years guarantee on the product

commencing on the date of purchase. Within

the guarantee period we will eliminate, free of

charge, any defects in the appliance resulting

from faults in materials or workmanship,

either by repairing or replacing the complete

appliance as we may choose. This guarantee

extends to every country where this appliance

is supplied by Braun or its appointed

distributor.

This guarantee does not cover: damage due to

improper use, normal wear or use, especially

regarding the brushheads, as well as defects

that have a negligible effect on the value or

operation of the appliance. The guarantee

becomes void if repairs are undertaken by

unauthorised persons and if original Braun

parts are not used.

To obtain service within the guarantee period,

hand in or send the complete appliance with

your sales receipt to an authorised Oral-B

Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights

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