SAVE THESE INSTRUCTIONS

- Never use while drowsy.
- This manual. Only use attachments recommended by the manufacturer.
- Use this product only for its intended use as described in used by, on, or near children or disabled persons.
- Close supervision is necessary when this product is used by, on, or near children or disabled persons.
- Keep the cord away from heated surfaces.
- the cord by grasping the plug – NOT THE CORD.
- DO NOT excessively bend or pull the cord. Unplug Refer to an Authorized Service Center for repair.
- or damaged, or if it has been immersed in water. if it is not working properly, if it has been dropped Never operate the unit if it has a damaged plug, other liquid.
- DO NOT submerge the unit in water or any 60 Hz outlet. Never force the plug into the outlet.
- The unit is designed for use in a 120 volt AC, Authorized Service Center after it has been in water. unplug immediately. Do not reach into the water under any circumstances. Have the unit checked by an
- DO NOT reach for a unit that has fallen into water.
- DO NOT pull into a tub, sink, or toilet.
- DO NOT place or store the unit where it can fall or ALWAYS unplug this product immediately after using.
- DO NOT use the unit while bathing.

DANGER
To reduce the risk of electrocution:

WARNING
To reduce the risk of burns, electrocution, fire, or injury:

SAVE THESE INSTRUCTIONS


- If you have had oral or periodontal surgery in the previous twelve months, contact your dentist prior to use.
- When you start any new oral hygiene program, your dentist if excessive bleeding occurs or if bleeding continues after two weeks of use.
- Be sure that only the water jet nozzle touches your teeth and gums. Damage to teeth and gums may occur if any other part of the unit touches them.
- This product is designed for oral cleaning only. DO NOT use the unit for any other purpose.
- As with any other new oral hygiene program, individuals with heart problems should check with their physician prior to using the Interplak® Dental Water Jet, and should always exercise particular care with any oral hygiene device.
- The Interplak® Dental Water Jet is a personal care unit and is intended for household use only.

MEDICAL WARNING
For your safety:

OPERATING INSTRUCTIONS

Cleaning and Care

If, in the course of use, the valve (8) on the water tank begins to leak, this can be taken care of by cleaning the valve seat. To do this, push out the valve body using the nozzle tip, clean the valve seat and replace the valve.

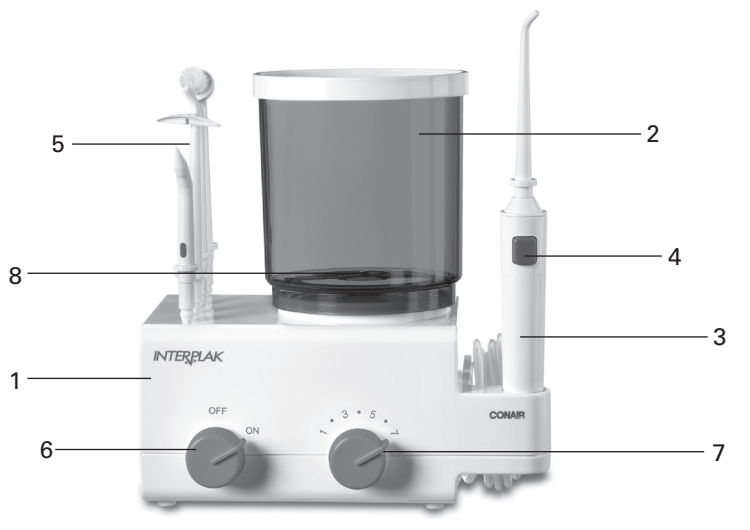- 1 Unit
- 2 Reservoir
- 3 Tip support
- 4 Finger tip on/off control
- 5 Color coded attachments
- 6 Power switch
- 7 Rotary flow control
- 8 Valve

Service Centers:
Conair Corporation
Service Department
150 Milford Road
East Windsor, NJ 08520
Glendale, AZ 85307

Conair Corporation
Service Department
7475 N. Glen Harbor Blvd.
East Windsor, NJ 08520

Limited 2-Year Warranty
Conair will repair or replace (at our option) your unit free of charge for 24 months from the date of purchase if the appliance is defective in workmanship or materials. To obtain service under this warranty, return the defective product to the service center listed below that is nearest you, together with your purchase receipt and \$3.00 for postage and handling. California residents need only provide proof of purchase and should call 1-800-334-4031 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 24 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 24 MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

WARNING
Changes or modifications to this unit not expressly approved by Conair could void the user's authority to operate the equipment and void limited warranty. Servicing should always be referred to an Authorized Service Center.

Limited 2-Year Warranty
Conair will repair or replace (at our option) your unit free of charge for 24 months from the date of purchase if the appliance is defective in workmanship or materials. To obtain service under this warranty, return the defective product to the service center listed below that is nearest you, together with your purchase receipt and \$3.00 for postage and handling. California residents need only provide proof of purchase and should call 1-800-334-4031 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 24 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 24 MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

Instructions For Use

Remove the water reservoir (2) and fill with lukewarm water. Reinstall the reservoir firmly onto the unit (1). Prior to initial use, the unit should be rinsed by letting a full reservoir of water pass through it. If the unit is to be used with a mouthwash solution, it should be rinsed with clean water after each use to remove any traces of mouthwash. If the unit has not been used for a long period of time, for reasons of hygiene it is advisable to follow the directions for initial use, letting a full reservoir of water pass through. To remove trapped air upon initial operation, turn the rotary flow control (7) all the way to the right and then to the left, three times.

Prior to filling the water reservoir, make sure water tip is firmly placed on and in the stand. The Water Jet is supplied with 5 attachments. Remove the tip support (3) and install selected tip:

- Standard Tip (2pcs):** removes trapped particles from between teeth.
- Subgingival Tip (1pc):** uses a soft rubber tip to massage the gums and safely deliver solution into the gingival pockets below the gumline.
- Brush Head Tip (1pc):** loosens debris and larger food particles before flushing them away – ideal for all dental work, including fixed and removable appliances.
- Tongue Cleaner (1pc):** place the tongue cleaner at farthest end of tongue and move forward. Repeat to cover entire tongue. After switching on (6) the unit, the rotary flow control (7) can be used to regulate the water pressure. You should begin with the lowest setting so that your gums can become accustomed to the water jet.


Hold the fingertip on/off control (water-stop button) (4) down and, after switching on the unit, place the tip in your mouth and lean over the sink. When you release the water-stop button, the cleaning and massage process is begun. The direction of the water jet in your mouth can be controlled by twisting the tip to the desired position on the handle. Keep your mouth slightly open and hold your elbow at mouth level so that the water can run back into the sink. Direct the water jet toward your teeth, especially between your teeth, in order to rinse away stubborn remaining food.

This also massages your gums. Do not point the water jet directly toward your gums; point between your teeth to prevent damage to delicate gum areas.

Welcome to a New Level of At-Home Dental Care

The Interplak® Water Jet enhances oral hygiene and significantly contributes to healthy gums when used as part of a responsible oral health management program. It is particularly effective at cleaning around orthodontic appliances, fixed and removable bridge work, implants and splints, and it provides valuable gum massage with every use.

Remember: Oral irrigation is a supplement to, not a substitute for, your oral hygiene routine. Ask your dental health care provider to advise you on how the Water Jet can most effectively meet your specific oral needs.

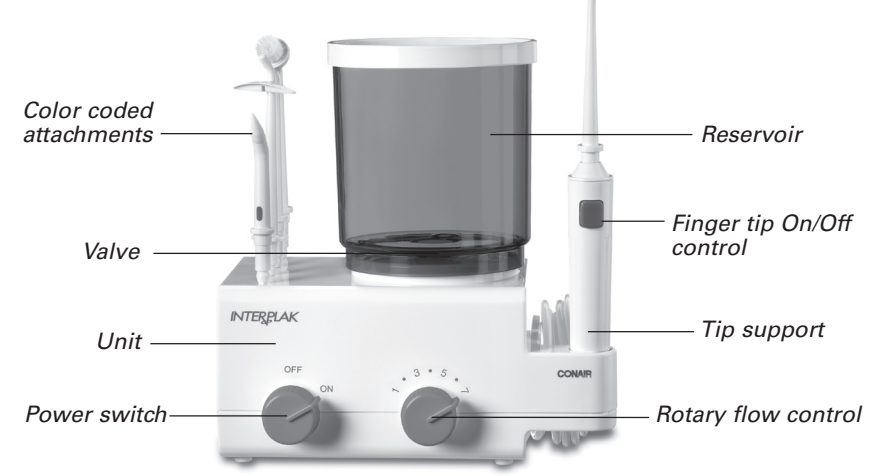- Color coded attachments
- Reservoir
- Finger tip On/Off control
- Tip support
- Rotary flow control
- Power switch
- Unit
- Valve

DENTAL WATER JET

Instruction Manual • Model WJB

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

Visit us at our website: www.interplak.com
Questions or Comments? 1-800-334-4031

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