

Health o meter®

P r o f e s s i o n a l

www.healthometermedical.com

800KL Instruction Manual

Please read this instruction manual carefully before use.

1. There is a plastic film with a red arrow under the Lithium battery cover to prevent battery drain in shipping. Please remove the strip before continuing with these instructions. Also, some models have a static cling label on the lens to prevent scratching. Please remove before use.
2. Your scale has been set at the factory to read weight in pounds. You may select kilograms or stone units of measurement by pressing the "kg/st/lb" button on the bottom of the scale.
3. Place scale on flat, hard surface. Carpeted or uneven floors may affect accuracy.
4. Step on the scale (disregarding the reading) and step off. This process insures the scale is set and is now ready for use. Do this whenever the scale is moved.
At all other times step straight onto the scale.
5. Step onto the scale platform and remain still. The display will show "0.0", then count up to your weight. The display will flash twice when the weight is locked. Your weight is now displayed.
6. If you step off of the scale, the display will turn off automatically within 5 seconds. If you continue to stand on the scale, your weight will be displayed for 10 seconds and then the scale will turn off automatically.

Precautions


1. Clean exposed parts with a soft, slightly damp cloth. To remove stains, use a mild soap. Never use detergents, excess water, treated cloths, harsh cleaning agents, or sprays. Do not attempt to lubricate or open scale casing.
2. Although your scale is designed to be maintenance-free, it contains sensitive electronic parts. Avoid rough treatment. Do not drop, kick, or jump on it. Treat it with care to ensure the best performance.
3. This scale is designed to read weights between 25 and 400 lbs / 12 and 180 kgs / 2 and 28 st 8 lbs. To prevent damage, do not place more than 400 lbs / 180 kgs / 28 st 8 lbs on it.
4. Do not store scale where you store cleaning chemicals. The vapors from some household products may affect the electronic components of your scale. Do not store scale on its side.
5. This scale is an extremely sensitive weighing device. In order to prevent run down of the battery, do not store anything on the scale.
6. Do not dispose of batteries in fire. Batteries may explode or leak. Remove the battery if the scale will not be used for a long period of time.

Battery Replacement

1. Use one 3 volt CR2032 lithium battery (as installed).
2. Open the battery compartment on the bottom of the scale by pushing the tab and lifting up the cover.
3. Remove the old battery from the compartment.
4. Place the new battery with the "+" sign up into the battery compartment.
5. Push gently down on the compartment cover until it clicks shut.

Low Battery

Your scale is equipped with a low battery indicator. It will display "LO" or readings will grow dim or become irregular if the battery needs replacement.

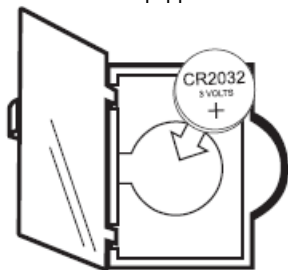

Trouble Shooting

1. Place your feet in the same position every time you weigh yourself. Stand still during measurement.
2. Weigh yourself without clothing whenever possible. Clothes can vary in weight.
3. Move the scale to a location where it won't be bumped. If your scale is bumped while being activated or while in use, an inaccurate reading may occur.
4. For the most consistent readings, weigh yourself at the same time of day. Allow for weight fluctuations due to food or drink in your system.
5. The maximum capacity of this scale is 400 lbs / 180 kgs / 28 st 8 lbs.. Should "Err" appear on the display during a weighing, this indicates that the maximum weighing capacity of the scale has been exceeded. Remove the weight immediately; otherwise, permanent damage may occur.

LIMITED WARRANTY

What does the warranty cover?

Pelstar LLC scales are warranted from date of purchase against defects of materials or in workmanship for a period of one (1) years. If scale fails to function properly, return the product, freight prepaid and properly packed to Pelstar. If manufacturer determines that a defect of material or in workmanship exists, customers' sole remedy will be repair or replacement of scale at no charge. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. All parts including repaired and replaced parts are covered only for the original warranty period.

Who is Covered?

The original purchaser of the product must have proof of purchase to receive warranty service. Pelstar dealers or retail stores selling Pelstar products do not have the right to alter, or modify or any way change the terms and conditions of this warranty.

What is Excluded?

Your warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instruction, abuse including tampering, damage in transit, or unauthorized repair or alternations. Further, the warranty does not cover Acts of God, such as fire, flood, hurricanes, and tornadoes. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country, state to state, province to province or jurisdiction to jurisdiction.

Warranty Service

Make sure you keep your sales receipt or document showing proof of purchase. Call 1 (800) 638-3722 or 1 (708) 598-9100 to receive a return authorization number. Attach proof of purchase to your defective product along with your name, address, daytime telephone number and description of the problem. Carefully package the product and send with shipping and insurance prepaid to: **Pelstar LLC Attention R/A# _____ Repair Department 7400 W. 100th Place Bridgeview, IL 60455** If your scale is not covered by warranty, or has been damaged, an estimate of repair costs, or replacement costs will be provided to you for approval prior to servicing or replacing.

Not legal for trade

Health o meter®
P r o f e s s i o n a l

7400 West 100th Place, Bridgeview IL 60455 • 1-800-815-6615

www.healthometermedical.com

Health o meter® is a registered trademark of Sunbeam Products Inc., Boca Raton, FL 33431

Free Manuals Download Website

<http://myh66.com>

<http://usermanuals.us>

<http://www.somanuals.com>

<http://www.4manuals.cc>

<http://www.manual-lib.com>

<http://www.404manual.com>

<http://www.luxmanual.com>

<http://aubethermostatmanual.com>

Golf course search by state

<http://golfingnear.com>

Email search by domain

<http://emailbydomain.com>

Auto manuals search

<http://auto.somanuals.com>

TV manuals search

<http://tv.somanuals.com>