

INVACARE® TURN-Q™ LTM LATERAL TURNING MATTRESS



The Invacare® Turn-Q™ LTM is a powered lateral turning mattress which automatically turns the user up to 30° bilaterally at preset intervals of 30, 60, or 90 minutes. The LTM reduces interface pressure below capillary closure for the treatment of stage I through stage IV pressure ulcers.

STANDARD FEATURES

- Automatic turning and static flotation modes
- Bilaterally turns patient up to 30° (depends on patient weight and body mass)
- Turn intervals 30, 60, and 90 minutes
- Power-failure mode keeps the lower mattress inflated for up to 12 hours in the event of a power failure
- Easy-to-understand, caregiver-friendly membrane control panel



INVACARE® TURN-Q™ LTM LATERAL TURNING MATTRESS



USER-FRIENDLY CONTROLS

- Power button with LED indicator displays whether the unit is in standby or operating mode
- Auto-firm with LED indicator inflates the mattress to maximum firmness to facilitate patient handling
- Patient weight adjustable from 70 - 350 lbs.
- Turning intervals of 30, 60, and 90 minutes, and a non-turning mode
- Auto-firm off will manually cancel the auto-firm mode; auto-firm mode will automatically return to normal operation if left unattended for 30 minutes or more




PORTABILITY

- Portable, lightweight mattress-replacement system fits most standard hospital beds
- Control unit easily attaches to top of footboard
- CPR quick-deflate handle activates rapid deflation of the mattress in approximately 30 seconds

PRESSURE RELIEF

Average pressure per pressure point for 10 subjects

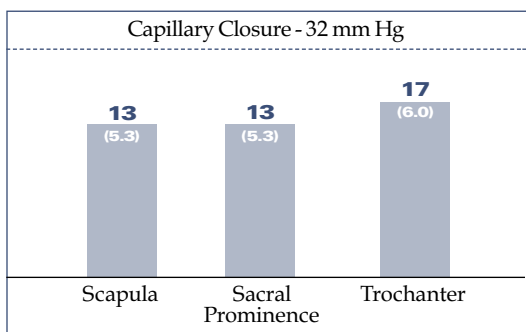

() Standard deviation

*Testing completed in an independent test lab

SPECIFICATIONS

POWER	115 VAC, 60 Hz, 15 watts
DIMENSIONS	
Standard mattress	35" W x 80" L x 11" H
Power unit	20.5" W x 10.8" L x 5.6" H
WEIGHT	
Mattress	29 lbs.
Power unit	23 lbs.
Maximum patient weight	350 lbs.
SAFETY	
Code approval	UL544
Fire	CAL 106 & 117
WARRANTY	
Power unit	1-year limited
Soft goods	6-month limited
MATERIAL	
Power unit	Cast, machined aluminum
Therapy cover	200-denier ripstop nylon with waterproof coating

Invacare Corporation www.invacare.com



Yes, you can.™

USA
One Invacare Way
Elyria, Ohio
44036-4028
(800) 333-6900

Canada
5970 Chedworth Way
Mississauga, Ontario
L5R 3T9 Canada
(800) 668-5324

Invacare, Yes, you can. and Turn-Q are trademarks of Invacare Corporation.
Specifications are subject to change without notification.

© 2000 Invacare Corporation Form No. 97-35 Rev. 05/00

Free Manuals Download Website

<http://myh66.com>

<http://usermanuals.us>

<http://www.somanuals.com>

<http://www.4manuals.cc>

<http://www.manual-lib.com>

<http://www.404manual.com>

<http://www.luxmanual.com>

<http://aubethermostatmanual.com>

Golf course search by state

<http://golfingnear.com>

Email search by domain

<http://emailbydomain.com>

Auto manuals search

<http://auto.somanuals.com>

TV manuals search

<http://tv.somanuals.com>