

The Invacare® Comet™ pediatric wheelchair's good looks are second only to its adjustability. The chair is designed to fit even the smallest children and grow right along with them, simply and inexpensively. A depth range of 7" to 13" is available on one frame. Or, choose the long frame for a growth range of 9" to 16". The chair's width can be changed just as easily, from 10" to 16" on either the short or long frame. Even the crossmembers can be moved to accommodate almost any seating system. For smaller children, the chair may be set up in a reversed configuration, with the larger wheels in front, and then switched to the conventional configuration when the child is able to reach the rear wheels — all with no additional parts required. These innovations have helped the Comet introduce a new word to the world of pediatric chairs: simplicity.

Highlights

- Swing-under footrest can be operated by the user or an attendant
- Reversible with no need for additional parts
- Available mini-tilt option offers 20° of tilt for children who do not need a full tilt-in-space system
- Wheels can be cambered
- 18" wheels available
- 16" wheels have true pediatric handrims
- Choice of short or long frame
- All rear wheel choices are no charge



The mini-tilt option is a simple and lightweight feature that provides up to 20° of tilt.



The Comet can be reversed with no additional parts required.



Available with swing-under footrest for easier transfers.

Specifications

Dimensions

Overall		
Width open	9" plus seat width*	
Front Seat-to-Floor Heights		
Caster & Wheel Configuration	Forward	Reverse
3" & 16"	10.75", 11.75"	10.75", 11.75"
3" & 18"	11.75", 12.75" ^{††} , 13.75"	11.75", 13.75"
3" & 20"	12.75" ^{††}	13.75", 14.5" ^{††} , 15" ^{†††}
3" & 22"	14.5", 15.5" ^{††}	14.5", 15.5" ^{††}
3" & 24"	15.5" ^{††} , 16.5" ^{†††}	16.5" ^{††}
5" & 16"	11.75"	16.5"
5" & 18"	11.75", 12.75" ^{††} , 13.75"	11.75", 12.75", 13.75"
5" & 20"	12.75", 13.75" ^{††}	13.75", 14.5" ^{††}
	14.5", 15.5" ^{††} , 16.5" ^{†††}	15.5" ^{††}
5" & 22"	14.5", 15.5"	14.5", 15.5" ^{††}
5" & 24"	15.5", 16.5" ^{†††}	16.5"
6" & 20"	16.5"	n/a
6" & 22"	16.5", 17.5"	n/a
6" & 24"	17.5" ^{††}	16.5", 17.5"

Knee-to-Heel** (range adjusts in 1" increments)

90° SWUNDER	4"-15.5" long frame 2"-8" short frame
90° FTBD	0"-15.5" long frame 0"-8" short frame
90° INDFPT	0"-16" long frame 0"-9" short frame
AT5046	12"-16" long frame n/a short frame

Arm Heights*** (range adjusts in 1" increments)

Cantilever	7"-12" add 1" for upholstered armpads
T-Arm	6"-10" add .5" for upholstered armpads 8"-14" add 1" for upholstered armpads

Cane Options

Adjustable-height	Sportster, 10° with push handles
Fixed-height	Straight cane 20"

*Measurement taken from outside of handrim to outside of handrim with 0° camber.

**Measurement taken from top of seat rail to back of footrest.

***Measurement taken from top of seat rail to top of armrest.

†Denotes 1" spacer.

††Denotes 2" spacer.

Seat-to-floor heights of 10.75", 11.75", 12.75", and 13.75" apply to 7"-13" adjustable seat depth only.

Seat-to-floor height increases by 1.75" with mini-tilt option.

Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

5970 Chedworth Way
Mississauga, Ontario
L5R 3T9 Canada
(800) 668-5324


Invacare is a registered trademark and Yes, you can. and Comet are trademarks of Invacare Corporation. Specifications are subject to change without notification.

© 2001

Invacare Corporation

Form No. 00-49 Rev. 04/01

Invacare® Comet™ Wheelchair

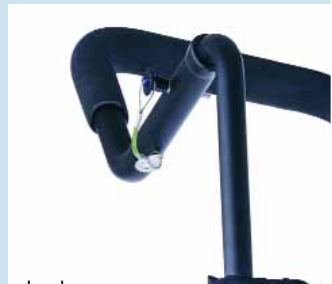

Free Manuals Download Website

<http://myh66.com>

<http://usermanuals.us>

<http://www.somanuals.com>

<http://www.4manuals.cc>

<http://www.manual-lib.com>

<http://www.404manual.com>

<http://www.luxmanual.com>

<http://aubethermostatmanual.com>

Golf course search by state

<http://golfingnear.com>

Email search by domain

<http://emailbydomain.com>

Auto manuals search

<http://auto.somanuals.com>

TV manuals search

<http://tv.somanuals.com>